

Email: info@biomatik.com https://www.biomatik.com

User Manual

Catalog Number: EKC31286

Product Name: Human microalbunminuria(MAU/ALB) ELISA kit

Detection Range: 0.078 μg/ml-5 μg/ml

Intended Use: For quantitative determination of human microalbunminuria

(MAU/ALB) concentrations in serum, plasma, urine.

Precautions: For research use only. Not for use in diagnostic procedures.

Manual Version: 202301V1

Storage:

Unopened kit	6 months when stored at 2 - 8°C.
Opened Kit	May be stored up to 1 month at 2 - 8°C. Keep it in sealed aluminum foil bag and avoid moisture.

The product manual may be updated as a result of continuous improvements.

Always refer to the hard copy manual included in the kit for your experiment.



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

Kit Components

Reagents	Quantity
Assay Plate (12 x 8 coated Microwells)	1(96 wells)
Standard (Freeze dried)	2
Biotin-conjugate (100 x concentrate)	1 x 60 µl
HRP-avdin (100 x concentrate)	1 x 120 µl
Biotin-conjugate Diluent	1 x 10 ml
HRP-avdin Diluent	1 x 20 ml
Sample Diluent	2 x 20 ml
Wash Buffer (25 x concentrate)	1 x 20 ml
TMB Substrate	1 x 10 ml
Stop Solution	1 x 10 ml
Adhesive Strip (For 96 wells)	4
Instruction Manual	1

Working Principle

This assay employs the competitive inhibition enzyme immunoassay technique. Antibody specific for MAU/ALB has been pre-coated onto a microplate. Standards and samples are pipetted into the wells with biotin-conjugated MAU/ALB. A competitive inhibition reaction is launched between MAU/ALB (Standards or samples) and biotin-conjugated MAU/ALB with the pre-coated MAU/ALB antibody. After washing, avidin conjugated Horseradish Peroxidase (HRP) is added to the wells. Following a wash to remove any unbound reagent, a substrate solution is added to the wells and color develops in opposite to the amount of MAU/ALB bound in the initial step. At last, measure the intensity of color after stopping color development.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Sensitivity

The minimum detectable dose of human MAU/ALB is less than 0.019

ug/ml. The sensitivity of this assay, or Lower Limit of Detection (LLD) was

defined as the lowest human MAU/ALB concentration that could be

differentiated from zero.

BIOMATIK

Specificity

This assay has high sensitivity and excellent specificity. No significant

cross-reactivity or interference between human MAU/ALB and analogues

was observed. Limited by current skills and knowledge, it is impossible

for us to complete the cross-reactivity detection between human

MAU/ALB and all the analogues, therefore, cross reaction may still exist.

Precision

Intra-assay Precision (Precision within an assay): CV%<8%

Three samples of known concentration were tested twenty times on one

plate to assess.

Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to

assess.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

BIOMATIK

Other Supplies Required

Microplate reader capable of measuring absorbance at 450 nm,

with the correction wavelength set at 600 nm or 630 nm.

An incubator which can provide stable incubation conditions up to

37°C±0.5°C.

Squirt bottle, manifold dispenser, or automated microplate washer.

Absorbent paper for blotting the microtiter plate.

100ml and 500ml graduated cylinders.

Deionized or distilled water.

Pipettes and pipette tips.

Test tubes for dilution.

Sample Collection & Storage

Serum Use a serum separator tube (SST) and allow samples to

clot for two hours at room temperature or overnight at 4°C before

centrifugation for 15 minutes at 1000 ×g. Remove serum and assay

immediately or aliquot and store samples at -20°C or -80°C. Avoid

repeated freeze-thaw cycles.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

BIOMATIK

Plasma Collect plasma using EDTA, or heparin as

anticoagulant. Centrifuge for 15 minutes at 1000 ×g at 2-8°C within 30

minutes of collection. Assay immediately or aliquot and store samples at -

20°C or -80°C. Avoid repeated freeze-thaw cycles.

Urine Use a sterile container to collect urine samples. Remove any

particulates by centrifugation for 15 minutes at 1000xg, 2 - 8°C and assay

immediately or aliquot and store samples at -20°C or -80°C. Avoid

repeated freeze-thaw cycles. Centrifuge again before assaying to remove

any additional precipitates that may appear after storage.

Sample Preparation

Recommend to dilute the serum or plasma samples 50000-fold before

test. The suggested 50000-fold dilution can be achieved by adding 2ul

sample to 398µl of normal saline. Complete the 50000-fold dilution by

adding 2µl of this solution to 498µl of Sample Diluent. The recommended

dilution factor is for reference only. The optimal dilution factor should be

determined by users according to their particular experiments.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Recommend to dilute the urine samples with Sample Diluent(1:40) before

test. The suggested 40-fold dilution can be achieved by adding 6µl

sample to 234ul of Sample Diluent. The recommended dilution factor is

for reference only. The optimal dilution factor should be determined by

users according to their particular experiments.

Notes:

1. Biomatik is only responsible for the kit itself, not for the samples

consumed during the assay. The user need to calculate the possible

amount of the samples to be used in the whole test. Please reserve

sufficient samples in advance.

BIOMATIK

2. Samples to be used within 5 days may be stored at 2-8°C, otherwise,

samples must be stored at -20°C (≤1month) or -80°C (≤2month) to avoid

contamination and loss of bioactivity.

3. Grossly hemolyzed samples are not suitable to use.

4. It would be necessary to run a preliminary experiment for validation, if

the samples are not indicated in the manual.

5. Please predict the concentration before assaying. If results were not

within the range of the standard curve, users would need to estimate the

optimal sample dilutions for their particular experiments.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

BIOMATIK

chemicals

6. Tissue or cell extraction samples prepared by chemical lysis buffer may cause unexpected ELISA results due to the impacts of certain

7. Considering the possibility of mismatch between antigen from other

resources and antibody in our kits (e.g., antibody targets conformational

epitope rather than linear epitope), some native or recombinant proteins

from other manufacturers may not be recognized by our products.

8. Due to factors including cell viability, cell number and sampling time.

samples from cell culture supernatant may not be detected by the kit.

9. Recommend to use fresh samples for the test. If you store samples for long time, protein degradation and denaturalization may occur in those

samples and finally lead to wrong results.

Reagents Preparation

Kindly use graduated containers to prepare the reagent.

Bring all reagents to room temperature (18-25°C) before use for 30

mins

Prepare fresh standard for each assay. Use within 4 hours and

discard after use

Making serial dilution in the wells directly is not permitted.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Please carefully reconstitute Standards according to the instruction,

and avoid foaming and mix gently until the crystals have completely

dissolved. To minimize imprecision caused by pipetting, use small

volumes and ensure that pipettors are calibrated. It is recommended to

suck more than 10_{ul} for once pipetting.

BIOMATIK

Distilled water is recommended to be used to make the preparation

for reagents. Contaminated water or container for reagent preparation will

influence the detection result

Biotin-conjugate (1x) - Centrifuge the vial before opening.

Biotin-conjugate requires a 100-fold dilution. A suggested 100-fold dilution

is 10 µl of Biotin-conjugate + 990 µl of Biotin-conjugate Diluent.

HRP-avidin (1x) - Centrifuge the vial before opening.

HRP-avidin requires a 100-fold dilution. A suggested 100-fold dilution is 10

μl of HRP-avidin + 990 μl of HRP-avidin Diluent.

Wash Buffer (1x)- If crystals have formed in the concentrate, warm them up

to room temperature and mix them gently until they get completely dissolved.

Dilute 20 ml of Wash Buffer Concentrate (25 x) into deionized or distilled

water to prepare 500 ml of Wash Buffer (1 x).




Email: info@biomatik.com

Standard

Centrifuge the standard vial at 6000-10000rpm for 30s.

Reconstitute the **Standard** with 1.0 ml of **Sample Diluent**. Do not substitute other diluents. This reconstitution produces a stock solution of 5 μ g/ml. Mix the standard to ensure complete reconstitution and allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions.

Pipette 150 μ l **Sample Diluent** into each tube (S0-S6). Use the stock solution to produce a 2-fold dilution series (below). Mix each tube thoroughly before the next transfer. The undiluted Standard serves as the high standard (5 μ g/ml). **Sample Diluent** serves as the zero standard (0 μ g/ml).

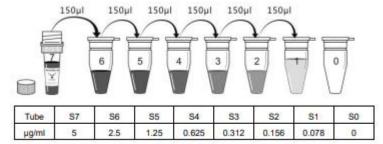

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Key Notes

Do not mix or substitute reagents with those from other lots or

sources

If samples generate higher values than the highest standard, dilute

the samples with Sample Diluent and repeat the assay.

Any variation in Sample Diluent, operator, pipetting technique,

washing technique, incubation time or temperature, and kit age can

cause variation in binding.

BIOMATIK

This assay is designed to eliminate interference by soluble

receptors, binding proteins, and other factors present in biological

samples. Until all factors have been tested in the Immunoassay, the

possibility of interference cannot be excluded.

Precautions

The Stop Solution provided with this kit is an acid solution. Wear eye, hand, face,

and clothing protection when using this material.

Assav Procedures

Bring all reagents and samples to room temperature before use.

Centrifuge the sample again after thawing before the assay. It is

recommended to assay all samples and standards in duplicate.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

BIOMATIK

1. Prepare all reagents, working standards, and samples as directed in

the previous sections.

2. Determine the number of wells to be used and put remaining wells and the

desiccant back into the pouch and seal the ziploc, store unused wells at 4°C.

3. Set a **Blank** well without any solution.

4. Add 50µl Standard or Sample per well.

5. Add 50µl Biotin-conjugate(1x) to each well (Not to the Blank well). Seal

plate with adhesive strip. Mix well and then incubate at 37°C for 60 mins.

(Biotin-conjugate (1x) may appear cloudy. Warm up to room temperature

and mix gently until solution appears uniform.)

6. Aspirate each well and wash, repeat the process two times for a total of

three washes. Wash by filling each well with Wash Buffer (200µl) using a

squirt bottle, multi-channel pipette, manifold dispenser, or autowasher, and

let it stand for 2 mins. Complete removal of liquid at each step is essential to

good performance. After the last wash, remove any remaining Wash Buffer

by aspirating or decanting. Invert the plate and blot it against clean paper

towels.

7. Add 100ul **HRP-avidin(1x)** to each well (not to Blank well). Seal plate with

a new adhesive strip. Mix well and then incubate at 37°C for 60 mins.

8. Repeat the aspiration/wash process for five times as in step 6.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

9. Add 90µl TMBSubstrate to each well and mix well. Cover the microtiter

plate adhesive strip. Incubate at 37°C for 20 mins. Protect from light.

10. Add 50ul Stop Solution to each well, gently tap the plate to ensure

thorough mixing.

BIOMATIK

11. Determine the optical density of each well within 5 mins, using a

microplate reader set to 450 nm. If wavelength correction is available, set to

540 nm or 570 nm. Subtract readings at 540 nm or 570 nm from the readings

at 450 nm. This subtraction will correct for optical imperfections in the plate.

Readings made directly at 450 nm without correction may be higher and less

accurate.

*Samples may require dilution. Please refer to Sample Preparation.

Notes:

The final experimental results will be closely related to validity of 1.

products, operation skills of end users and the experimental

environments

2 Samples or reagents addition: Please use the freshly prepared

Standard. Please carefully add samples to wells and mix gently to avoid

foaming. Do not touch the well wall as possible. For each step in the

procedure, total dispensing time for addition of reagents or samples to

the assay plate should not exceed 10 mins. This will ensure equal



Email: info@biomatik.com https://www.biomatik.com



elapsed time for each pipetting step, without interruption. Duplication of all standards and specimens, although not required, is recommended. To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.

- 3. Incubation: To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary. Do not allow wells to sit uncovered for extended periods between incubation steps. Once reagents have been added to the well strips, DO NOT let the strips DRY at any time during the assay. Incubation time and temperature must be observed.
- Washing: The wash procedure is critical. Complete removal of 4. liquid at each step is essential to good performance. After the last wash. remove any remaining Wash Solution by aspirating or decanting and remove any drop of water and fingerprint on the bottom of the plate. Insufficient washing will result in poor precision and falsely elevated absorbance reading. When using an automated plate washer, adding a 30 second soak period following the addition of wash buffer, and/or rotating the plate 180 degrees between wash steps may improve assay precision.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

BIOMATIK

5. Controlling of reaction time: Observe the change of color after

adding TMB Substrate (e.g. observation once every 10 mins), TMB

Substrate should change from colorless or light blue to gradations of blue.

If the color is too deep, add Stop Solution in advance to avoid excessively

strong reaction which will result in inaccurate absorbance reading.

TMB Substrate is easily contaminated. TMB Substrate should 6.

remain colorless or light blue until added to the plate. Please protect it

from light.

7. Stop Solution should be added to the plate in the same order as

the TMB Substrate. The color developed in the wells will turn from blue to

yellow upon addition of the Stop Solution. Wells that are green in color

indicate that Stop Solution has not mixed thoroughly with TMB Substrate.

Calculation of Results

Using the professional soft "Curve Expert" to make a standard

curve is recommended, which can be downloaded from our web.

Average the duplicate readings for each standard and sample and

subtract the average optical density of Blank.





Email: info@biomatik.com https://www.biomatik.com

BIOMATIK

Create a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the x-axis against the concentration on the y-axis and draw a best fit curve through the points on the graph. The data may be linearized by plotting the log of the MAU/ALB concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data.

If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.