

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

User Manual

Catalog Number: EKC31578

Product Name: Bovine Caspase 7 (Casp-7) ELISA Kit

Detection Range: 0.625 ng/mL- 40 ng/mL

Intended Use: For quantitative determination of bovine caspase 7 (Casp-7)

concentrations in serum, plasma, tissue homogenates.

Precautions: For research use only. Not for use in diagnostic procedures.

Manual Version: 202301V1

Storage:

Unopened kit	6 months when stored at 2 - 8°C.
Opened Kit	May be stored up to 1 month at 2 - 8°C. Keep it in sealed aluminum foil bag and avoid moisture.

The product manual may be updated as a result of continuous improvements.

Always refer to the hard copy manual included in the kit for your experiment.



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140. (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

Kit Components

Reagents	Quantity
Assay plate (12 x 8 coated Microwells)	1(96 wells)
Standard (Freeze dried)	2
Biotin-antibody (100 x concentrate)	1 x 120 µl
HRP-avidin (100 x concentrate)	1 x 120 µl
Biotin-antibody Diluent	1 x 15 ml
HRP-avidin Diluent	1 x 15 ml
Sample Diluent	1 x 50 ml
Wash Buffer (25 x concentrate)	1 x 20 ml
TMB Substrate	1 x 10 ml
Stop Solution	1 x 10 ml
Adhesive Strip (For 96 wells)	4
Instruction manual	1

Working Principle

This assay adopts quantitative sandwich enzyme immunoassay technique. Antibody specific to Casp-7 has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any Casp-7 present is bound by the immobilized antibody. After removing any unbound substances, a biotin-conjugated antibody specific for Casp-7 is added to the wells. After washing, avidin conjugated Horseradish Peroxidase (HRP) is added to the wells. Following a wash to remove any unbound avidin-enzyme reagent, a substrate solution is added to the wells and color develops in proportion to the amount of Casp-7 bound in the initial step. At last, measure the intensity of color after stopping color development.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Sensitivity

BIOMATIK

The minimum detectable dose of bovine Casp-7 is less than 0.156 ng/mL.

The sensitivity of this assay, or Lower Limit of Detection (LLD) was

defined as the lowest protein concentration that could be differentiated

from zero. It was determined by the mean O.D value of 20 replicates of

the zero standard added by their three standard deviations.

Specificity

This assay has high sensitivity and excellent specificity. No significant

cross-reactivity or interference between bovine Casp-7 and analogues

was observed. Limited by current skills and knowledge, it is impossible

for us to complete the cross-reactivity detection between bovine Casp-7

and all the analogues, therefore, cross reaction may still exist.

Precision

Intra-assay Precision (Precision within an assay): CV%<8%

Three samples of known concentration were tested twenty times on one

plate to assess.

Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to

assess.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Other Supplies Required

BIOMATIK

Microplate reader capable of measuring absorbance at 450 nm,

with the correction wavelength set at 600 nm or 630 nm.

An incubator which can provide stable incubation conditions up to

37°C±0.5°C.

Squirt bottle, manifold dispenser, or automated microplate washer.

Absorbent paper for blotting the microtiter plate.

100ml and 500ml graduated cylinders.

Deionized or distilled water.

Pipettes and pipette tips.

Test tubes for dilution.

Sample Collection & Storage

Serum Use a serum separator tube (SST) and allow samples to

clot for two hours at room temperature or overnight at 4°C before

centrifugation at 1000 ×g for 15 mins. Remove serum and assay

immediately or aliquot and store samples at -20°C or -80°C. Avoid

repeated freeze-thaw cycles.



- Plasma Collect plasma using EDTA or heparin as an anticoagulant. Within 30 mins after collecting samples, centrifuge samples at 1000 x g, 2 - 8°C, for 15 mins. Assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles.
- **Tissue Homogenates** 100 mg tissue was rinsed with 1X PBS, homogenized in 1 ml of 1X PBS and stored overnight at -20°C. After two freeze-thaw cycles were performed to break the cell membranes, the homogenates were centrifuged for 5 minutes at 5000 x g. 2 - 8°C. The supernate was removed and assayed immediately. Alternatively, aliquot and store samples at -20°C or -80°C. Centrifuge the sample again after thawing before the assay. Avoid repeated freeze-thaw cycles.

Notes:

- 1. Biomatik is only responsible for the kit itself, not for the samples consumed during the assay. The user need to calculate the possible amount of the samples to be used in the whole test. Please reserve sufficient samples in advance.
- 2. Samples to be used within 5 days may be stored at 2-8°C, otherwise, samples must be stored at -20°C (≤1month) or -80°C (≤2month) to avoid contamination and loss of bioactivity.
- 3. Grossly hemolyzed samples are not suitable to use.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com
https://www.biomatik.com

BIOMATIK

4. It would be necessary to run a preliminary experiment for validation, if the samples are not indicated in the manual.

5. Please predict the concentration before assaying. If results were not within the range of the standard curve, users would need to estimate the optimal sample dilutions for their particular experiments.

Tissue or cell extraction samples prepared by chemical lysis buffer may cause unexpected ELISA results due to the impacts of certain chemicals.

7. Considering the possibility of mismatch between antigen from other resources and antibody in our kits (e.g., antibody targets conformational epitope rather than linear epitope), some native or recombinant proteins from other manufacturers may not be recognized by our products.

8. Due to factors including cell viability, cell number and sampling time, samples from cell culture supernatant may not be detected by the kit.

Recommend to use fresh samples for the test. If you store samples for long time, protein degradation and denaturalization may occur in those samples and finally lead to wrong results.



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

Reagents Preparation

- Kindly use graduated containers to prepare the reagent.
- Bring all reagents to room temperature (18-25°C) before use for 30

mins.

- Prepare fresh standard for each assay. Use within 4 hours and discard after use.
- Making serial dilution in the wells directly is not permitted.
- Please carefully reconstitute Standards according to the instruction, and avoid foaming and mix gently until the crystals have completely dissolved.
 To minimize imprecision caused by pipetting, use small volumes and ensure that pipettors are calibrated. It is recommended to suck more than 10 µl for once pipetting.
- Distilled water is recommended to be used to make the preparation for reagents or samples. Contaminated water or container for reagent preparation will influence the detection result.

Biotin-antibody (1x) - Centrifuge the vial before opening.

Biotin-antibody requires a 100 -fold dilution. A suggested 100 -fold dilution is 10 µl of Biotin-antibody + 990 µl of Biotin-antibody Diluent.

HRP-avidin (1x) - Centrifuge the vial before opening.
HRP-avidin requires a 100 -fold dilution. A suggested 100 -fold dilution is
10 μl of HRP-avidin + 990 μl of HRP-avidin Diluent.

Wash Buffer (1x) - If crystals have formed in the concentrate, warm up to room temperature and mix gently until the crystals have completely dissolved. Dilute 20 ml of Wash Buffer Concentrate (25 x) into deionized or distilled water to prepare 500 ml of Wash Buffer (1 x).



Tel: (519) 489-7195, (800) 836-8089


Fax: (519) 231-0140, (877) 221-3515 Email: info@biomatik.com https://www.biomatik.com

Standard

Centrifuge the standard vial at 6000 - 10000rpm for 30s.

Reconstitute the **Standard** with 1.0 ml of **Sample Diluent**. Do not substitute other diluents. This reconstitution produces a stock solution of 40 ng/ml. Mix the standard to ensure complete reconstitution and allow the standard to sit for a minimum of 15 mins with gentle agitation prior to making dilutions.

Pipette 250 μl **Sample Diluent** into each tube (S0-S6). Use the stock solution to produce a 2 -fold dilution series (below). Mix each tube thoroughly before the next transfer. The undiluted Standard serves as the high standard (40 ng/ml). **Sample Diluent** serves as the zero standard (0 ng/ml).

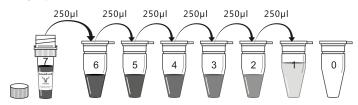

Tube	S7	S6	S5	S4	S3	S2	S1	S0
ng/ml	40	20	10	5	2.5	1.25	0.625	0

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Key Notes

Do not mix or substitute reagents with those from other lots or

sources

If samples generate higher values than the highest standard, dilute

the samples with Sample Diluent and repeat the assay.

Any variation in Sample Diluent, operator, pipetting technique,

washing technique, incubation time or temperature, and kit age can

cause variation in binding.

BIOMATIK

This assay is designed to eliminate interference by soluble

receptors, binding proteins, and other factors present in biological

samples. Until all factors have been tested in the Immunoassay, the

possibility of interference cannot be excluded.

Precautions

The Stop Solution provided with this kit is an acid solution. Wear eye, hand, face,

and clothing protection when using this material.

Assav Procedures

Bring all reagents and samples to room temperature before use.

Centrifuge the sample again after thawing before the assay. It is

recommended to assay all samples and standards in duplicate.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

1. Prepare all reagents and samples as directed in the previous sections.

2. Refer to the Assay Layout Sheet to determine the number of wells to be

used and put any remaining wells and the desiccant back into the pouch and

seal the ziploc, store unused wells at 4°C.

® BIOMATIK

3. Add 100µl **Standard** or **Sample** per well. Seal plate with the adhesive strip.

Incubate at 37°C for 2 hours. A plate layout is provided to record standards

and samples assayed.

4. Remove the liquid of each well, don't wash.

5. Add 100 µl Biotin-antibody (1x) to each well. Seal plate with a new

adhesive strip. Incubate at 37°C for 1 hour. (Biotin-antibody (1x) may

appear cloudy. Warm up to room temperature and mix gently until solution

appears uniform.)

6. Aspirate each well and wash, repeating the process two times for a total of

three washes. Wash by filling each well with Wash Buffer (200ul) using a

squirt bottle, multi-channel pipette, manifold dispenser, or autowasher, and

let it stand for 2 mins, complete removal of liquid at each step is essential to

good performance. After the last wash, remove any remaining Wash Buffer

by aspirating or decanting. Invert the plate and blot it against clean paper

towels.

7. Add 100ul **HRP-avidin (1x)** to each well. Seal plate with a new adhesive

strip and then incubate at 37°C for 60 mins.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

8. Repeat the aspiration/wash process for five times as in step 6.

9. Add 90ul TMB Substrate to each well, mix well, Incubate at 37°C for 15-

30 mins. Protect from light.

® BIOMATIK

10. Add 50µl Stop Solution to each well, gently tap the plate to ensure

thorough mixing.

11. Determine the optical density of each well within 5 mins, using a microplate

reader set to 450 nm.lf wavelength correction is available, set to 540 nm or 570

nm. Subtract readings at 540 nm or 570 nm from the readings at 450 nm. This

subtraction will correct for optical imperfections in the plate. Readings made

directly at 450 nm without correction may be higher and less accurate.

*Samples may require dilution. Please refer to Sample Preparation section.

Notes:

1. The final experimental results will be closely related to validity of

the products, operation skills of the end users and the experimental

environments

2 Samples or reagents addition: Please use the freshly prepared

Standard, Please carefully add samples to wells and mix gently to avoid

foaming. Do not touch the well wall as possible. For each step in the

procedure, total dispensing time for addition of reagents or samples to

the assay plate should not exceed 10 mins. This will ensure equal



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

BIOMATIK

elapsed time for each pipetting step, without interruption. Duplication of all standards and specimens, although not required, is recommended. To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions.

Also, use separate reservoirs for each reagent.

3. Incubation: To ensure accurate results, proper adhesion of plate

sealers during incubation steps is necessary. Do not allow wells to sit

uncovered for extended periods between incubation steps. Once

reagents have been added to the well strips, DO NOT let the strips DRY

at any time during the assay. Incubation time and temperature must be

observed.

Washing: The wash procedure is critical. Complete removal of 4.

liquid at each step is essential to good performance. After the last wash.

remove any remaining Wash Solution by aspirating or decanting and

remove any drop of water and fingerprint on the bottom of the plate.

Insufficient washing will result in poor precision and falsely elevated

absorbance reading. When using an automated plate washer, adding a

30 second soak period following the addition of wash buffer, and/or

rotating the plate 180 degrees between wash steps may improve assay

precision.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

BIOMATIK

5. Controlling of reaction time: Observe the change of color after

adding TMB Substrate (e.g. observation once every 10 mins), TMB

Substrate should change from colorless or light blue to gradations of blue.

If the color is too deep, add Stop Solution in advance to avoid excessively

strong reaction which will result in inaccurate absorbance reading.

TMB Substrate is easily contaminated. TMB Substrate should 6.

remain colorless or light blue until added to the plate. Please protect it

from light.

7. Stop Solution should be added to the plate in the same order as

the TMB Substrate. The color developed in the wells will turn from blue to

yellow upon addition of the Stop Solution. Wells that are green in color

indicate that Stop Solution has not mixed thoroughly with TMB Substrate.

Calculation of Results

Using the professional soft "Curve Expert" to make a standard

curve is recommended, which can be downloaded from our web.

Average the duplicate readings for each standard and sample and

subtract the average zero standard optical density.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Create a standard curve by reducing the data using computer software

capable of generating a four parameter logistic (4-PL) curve-fit. As an

alternative, construct a standard curve by plotting the mean absorbance

for each standard on the x-axis against the concentration on the y-axis

and draw a best fit curve through the points on the graph. The data may

be linearized by plotting the log of the Casp-7 concentrations versus the

log of the O.D. and the best fit line can be determined by regression

analysis. This procedure will produce an adequate but less precise fit of

the data.

If samples have been diluted, the concentration read from the standard

curve must be multiplied by the dilution factor.

BIOMATIK