

Biomatik Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515 Email: info@biomatik.com http://www.biomatik.com

CAU21580

Polyclonal Antibody to Hairy And Enhancer Of Split 1 (HES1)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)



Biomatik Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515 Email: info@biomatik.com http://www.biomatik.com

## [ PROPERTIES ]

Source: Polyclonal antibody preparation

Host: Rabbit

Purification: Antigen-specific affinity chromatography followed by Protein A affinity

chromatography

Traits: Liquid

Concentration: 1mg/mL

**UOM:** 100ul

Applications: WB; IHC; ICC; IP.

#### [IMMUNOGEN]

Immunogen: Recombinant HES1 (Val21~Asn251) expressed in E.coli

Accession No.: RPK289Hu01

### [APPLICATIONS]

Western blotting: 0.5-2µg/mL;1:500-2000

Immunohistochemistry: 5-20µg/mL;1:50-200

Immunocytochemistry: 5-20µg/mL;1:50-200

Optimal working dilutions must be determined by end user.

#### [FORMULATION]

**Form & Buffer:** Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

#### [STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 4°C for frequent use.

Aliquot and store at -20°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the



expiration date under appropriate storage condition.

### [IDENTIFICATION]

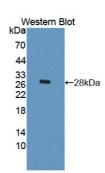


Figure. Western Blot; Sample: Recombinant HES1, Human.

# [<u>IMPORTANT NOTE</u>]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.

Biomatik Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515 Email: info@biomatik.com http://www.biomatik.com