

CAU25780 Mouse Polyclonal Antibody to Prominin 1 (PROM1) Organism Species: *Mus musculus* (Mouse) Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)

Biomatik - One Stop for Bio-Reagents, Custom Gene, Peptide, Protein and Antibody Services Tel: (800) 836-8089, (519) 489-7195 Fax: (519) 231-0140 www.biomatik.com



[PROPERTIES]

Source: Polyclonal antibody preparation Host: Rabbit Purification: Antigen-specific Affinity Chromatography. Traits: Liquid Concentration: 200µg/mL UOM: 100µg Applications: WB; IHC; ICC; IP.

[<u>IMMUNOGEN</u>]

Immunogen: Recombinant PROM1 (Ala180~Trp435) expressed in *E.coli*. Accession No.: RPB516Mu01

[APPLICATIONS]

Western blotting: 0.5-2µg/mL Immunohistochemistry: 5-20µg/mL Immunocytochemistry: 5-20µg/mL Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

[QUALITY CONTROL]

Content: The quality control contains recombinant PROM1 disposed in loading buffer.

 $\label{eq:Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.$

5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.



Loading Buffer: 100mM Tris(pH6.8), 1% SDS, 150mM NaCl, 50% glycerol, 0.02% BPB, 50mM DTT and 0.02% NaN₃.


[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.


Store at 4°C for frequent use.

Aliquot and store at -20°C for 12 months.


Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

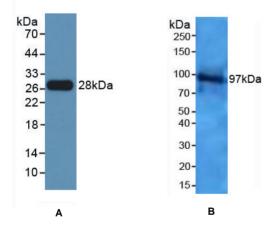

[IDENTIFICATION]

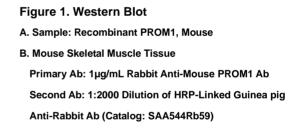



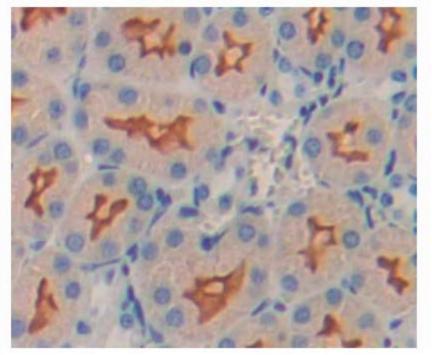

Used in DAB staining on fromalin fixed

paraffin- embedded kidney tissue