CAU27228 Polyclonal Antibody to Toll Like Receptor 4 (TLR4) Organism Species: Mus musculus (Mouse) Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)

[PROPERTIES]

Source: Polyclonal antibody preparation Host: Rabbit Purification: Antigen-specific Affinity Chromatography. Traits: Liquid Concentration: 200µg/mL UOM: 10µg Applications: WB; ICC; IHC-P; IHC-F; ELISA; IP; IF; FCM.

[IMMUNOGEN]

Immunogen: Recombinant TLR4 (Asn26~Glu269) expressed in *E.coli*. Accession No.: RPA753Mu01

[APPLICATIONS]

Western blotting: 0.5-2ug/ml Immunocytochemistry in formalin fixed cells: 5-20ug/ml Immunohistochemistry in formalin fixed frozen section: 5-20ug/ml Immunohistochemistry in paraffin section: 5-20ug/ml Enzyme-linked Immunosorbent Assay: 0.05-2ug/ml Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

[QUALITY CONTROL]

Content: The quality control contains recombinant TLR4 disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes. **Loading Buffer:** 100mM Tris(pH6.8), 1% SDS, 150mM NaCl, 50% glycerol, 0.02% BPB, 50mM DTT and 0.02% NaN₃.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 4°C for frequent use.

Aliquot and store at -20°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

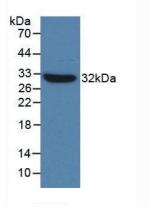


Figure 1. Western Blot Sample: Recombinant TLR4, Mouse

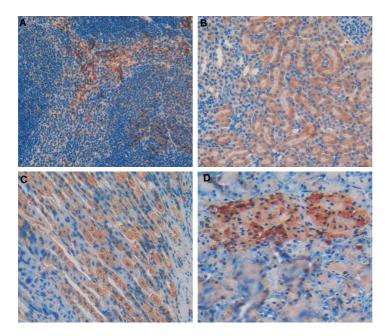


Figure 2. DAB staining on IHC-P

Samples:

- A. Mouse Spleen Tissue
- B. Mouse Kidney Tissue
- C. Mouse Stomach Tissue
- D. Mouse Pancreas Tissue