

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140. (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

User Manual

Catalog Number: EKC31504

Product Name: Sheep Albumin (Alb) ELISA Kit

Detection Range: 21.88 µg/ml - 350 µg/ml

Intended Use: For quantitative determination of sheep albumin (Alb)

concentrations in serum, plasma, cell culture supernates, urine, tissue

homogenates.

Precautions: For research use only. Not for use in diagnostic procedures.

Manual Version: 202301V1

Storage:

Unopened kit	6 months when stored at 2 - 8°C.
Opened Kit	May be stored up to 1 month at 2 - 8°C. Keep it in sealed
	aluminum foil bag and avoid moisture.

The product manual may be updated as a result of continuous improvements.

Always refer to the hard copy manual included in the kit for your experiment.





Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

Kit Components

Reagents	Quantity
Assay Plate (12 x 8 coated Microwells)	1(96 wells)
Standard	2 x 500 µl
HRP-conjugate (100 x concentrate)	1 x 60 µl
HRP-conjugate Diluent	1 x 10 ml
Sample Diluent	2 x 20 ml
Wash Buffer (25 x concentrate)	1 x 20 ml
TMB Substrate	1 x 10 ml
Stop Solution	1 x 10 ml
Adhesive Strip (For 96 wells)	4
Instruction Manual	1

Working Principle

This assay employs the competitive inhibition enzyme immunoassay technique. Antibody specific for Alb has been pre-coated onto a microplate. Standards or samples are added to the appropriate microtiter plate wells with Horseradish Peroxidase (HRP) conjugated antibody specific to Alb. The competitive inhibition reaction is launched between with pre-coated Alb and Alb in samples. A substrate solution is added to the wells and the color develops in opposite to the amount of Alb in the sample. At last, measure the intensity of color after stopping color development.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Sensitivity

BIOMATIK

The minimum detectable dose of sheep Alb is less than 5.5 µg/ml. The

sensitivity of this assay, or Lower Limit of Detection (LLD) was defined as

the lowest sheep Alb concentration that could be differentiated from zero.

Specificity

This assay has high sensitivity and excellent specificity. No significant

cross-reactivity or interference between sheep Alb and analogues was

observed. Limited by current skills and knowledge, it is impossible for us

to complete the cross-reactivity detection between sheep Alb and all the

analogues, therefore, cross reaction may still exist.

Precision

Intra-assay Precision (Precision within an assay): CV%<12%

Three samples of known concentration were tested twenty times on one

plate to assess.

Inter-assay Precision (Precision between assays): CV%<15%

Three samples of known concentration were tested in twenty assays to

assess.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Other Supplies Required

BIOMATIK

Microplate reader capable of measuring absorbance at 450 nm,

with the correction wavelength set at 600 nm or 630 nm.

An incubator which can provide stable incubation conditions up to

37°C±0.5°C.

Squirt bottle, manifold dispenser, or automated microplate washer.

Absorbent paper for blotting the microtiter plate.

100ml and 500ml graduated cylinders.

Deionized or distilled water.

Pipettes and pipette tips.

Test tubes for dilution.

Sample Collection & Storage

Serum Use a serum separator tube (SST) and allow samples to

clot for two hours at room temperature or overnight at 4°C before

centrifugation at 1000 ×g for 15 mins. Remove serum and assay

immediately or aliquot and store samples at -20°C or -80°C. Avoid

repeated freeze-thaw cycles.



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

Plasma Collect plasma using EDTA or heparin as an anticoagulant.
 Within 30 mins after collecting samples, centrifuge samples at 1000 x g, 2
 8°C, for 15 mins. Assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles.

- Cell Culture Supernates Remove particulates by centrifugation for 15 minutes at 1000 x g, 2 8°C and assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles.
- Urine Use a sterile container to collect urine samples. Remove
 any particulates by centrifugation for 15 minutes at 1000xg, 2 8°C and
 assay immediately or aliquot and store samples at -20°C or -80°C. Avoid
 repeated freeze-thaw cycles. Centrifuge again before assaying to remove
 any additional precipitates that may appear after storage.
- Tissue Homogenates 100mg tissue was rinsed with 1X PBS, homogenized in 1 ml of 1X PBS and stored overnight at -20°C. After two freeze-thaw cycles were performed to break the cell membranes, the homogenates were centrifuged for 5 minutes at 5000 x g, 2 8°C. The supernate was removed and assayed immediately. Alternatively, aliquot and store samples at -20°C or -80°C. Centrifuge the sample again after thawing before the assay. Avoid repeated freeze-thaw cycles.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Sample Preparation

BIOMATIK

Recommend to dilute the serum or plasma samples with Sample

Diluent(1:200) before test. The suggested 200-fold dilution can be

achieved by adding 5µl sample to 45µl of Sample Diluent. Complete the

200-fold dilution by adding15µl of this solution to 285µl of Sample Diluent.

The recommended dilution factor is for reference only. The optimal

dilution factor should be determined by users according to their particular

experiments.

Notes:

1. Biomatik is only responsible for the kit itself, not for the samples

consumed during the assay. The user need to calculate the possible

amount of the samples to be used in the whole test. Please reserve

sufficient samples in advance.

2. Samples to be used within 5 days may be stored at 2-8°C, otherwise,

samples must be stored at -20°C (≤1month) or -80°C (≤2month) to avoid

contamination and loss of bioactivity.

3. Grossly hemolyzed samples are not suitable to use.

4. It would be necessary to run a preliminary experiment for validation, if

the samples are not indicated in the manual.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

5. Please predict the concentration before assaying. If results were not

within the range of the standard curve, users would need to estimate the

optimal sample dilutions for their particular experiments.

BIOMATIK

6. Tissue or cell extraction samples prepared by chemical lysis buffer

may cause unexpected ELISA results due to the impacts of certain

chemicals.

7. Considering the possibility of mismatch between antigen from other

resources and antibody in our kits (e.g., antibody targets conformational

epitope rather than linear epitope), some native or recombinant proteins

from other manufacturers may not be recognized by our products.

8. Due to factors including cell viability, cell number and sampling time,

samples from cell culture supernatant may not be detected by the kit.

9. Recommend to use fresh samples for the test. If you store samples for

long time, protein degradation and denaturalization may occur in those

samples and finally lead to wrong results.

Reagents Preparation

Kindly use graduated containers to prepare the reagent.

Bring all reagents to room temperature (18-25°C) before use for 30

mins.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

BIOMATIK

Prepare fresh standard for each assay. Use within 4 hours and discard

after use

It is not allowed to make serial dilution in the wells directly.

Distilled Water is recommended. Contaminated water or container for

reagents preparation will affect the test result.

HRP-conjugate (1x) - Centrifuge the vial before opening.

HRP-conjugate requires a 100-fold dilution. A suggested 100-fold dilution is

10 µl of HRP-conjugate + 990 µl of HRP-conjugate Diluent.

Wash Buffer (1x)- If crystals have formed in the concentrate, warm them up

to room temperature and mix them gently until they get completely dissolved.

Dilute 20 ml of Wash Buffer Concentrate (25 x) into deionized or distilled

water to prepare 500 ml of Wash Buffer (1 x).

Standard

Centrifuge the standard vial at 6000-10000rpm for 30s.

Allow the stock standard solution of 350µg/ml to sit for a minimum of 30 mins to


warm to room temperature prior to making dilutions.



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

Pipette 150 μl **Sample Diluent** into each tube (S0-S4). Use the **undiluted Standard** (S5) solution to produce a 2-fold dilution series (below). Mix each tube thoroughly before the next transfer. **Sample Diluent** serves as the zero standard (0 μg/ml)).

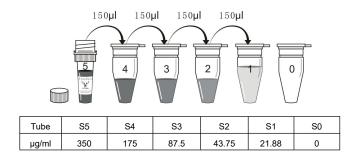

Key Notes

- Do not mix or substitute reagents with those from other lots or sources.
- If samples generate higher values than the highest standard, dilute the samples with Sample Diluent and repeat the assay.
- Any variation in Sample Diluent, operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

BIOMATIK https://www.biomatik.com

This assay is designed to eliminate interference by soluble

receptors, binding proteins, and other factors present in biological

samples. Until all factors have been tested in the Immunoassay, the

possibility of interference cannot be excluded.

Precautions

The Stop Solution provided with this kit is an acid solution. Wear eye, hand, face,

and clothing protection when using this material.

Assay Procedures

Bring all reagents and samples to room temperature before use.

Centrifuge the sample again after thawing before the assay. It is

recommended to assay all samples and standards in duplicate.

1. Prepare all reagents, working standards, and samples as directed in

the previous sections.

2. Determine the number of wells to be used and put remaining wells and the

desiccant back into the pouch and seal the ziploc, store unused wells at 4°C.

3. Set a **Blank** well without any solution.

4. Add 50µl Standard or Sample per well.

5. Add 50µl HRP-conjugate(1x) to each well (Not to the Blank well). Seal

plate with adhesive strip. Mix well and then incubate at 37°C for 40 mins. A

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

plate layout is provided to record standards and samples assayed.

6. Aspirate each well and wash, repeat the process four times for a total of

five washes. Wash by filling each well with Wash Buffer (200ul) using a

squirt bottle, multi-channel pipette, manifold dispenser, or autowasher, and

let it stand for 2 mins. Complete removal of liquid at each step is essential to

good performance. After the last wash, remove any remaining Wash Buffer

by aspirating or decanting. Invert the plate and blot it against clean paper

towels

7. Add 90µl TMBSubstrate to each well and mix well. Incubate at 37°C for

20 mins. Protect from light.

BIOMATIK

8. Add 50µl Stop Solution to each well, gently tap the plate to ensure

thorough mixing.

9. Determine the optical density of each well within 5 mins, using a

microplate reader set to 450 nm. If wavelength correction is available, set to

540 nm or 570 nm. Subtract readings at 540 nm or 570 nm from the readings

at 450 nm. This subtraction will correct for optical imperfections in the plate.

Readings made directly at 450 nm without correction may be higher and less

accurate.

*Samples may require dilution. Please refer to Sample Preparation.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Notes:

The final experimental results will be closely related to validity of 1

products, operation skills of end users and the experimental

environments.

BIOMATIK

2. Samples or reagents addition: Please use the freshly prepared

Standard. Please carefully add samples to wells and mix gently to avoid

foaming. Do not touch the well wall as possible. For each step in the

procedure, total dispensing time for addition of reagents or samples to

the assay plate should not exceed 10 mins. This will ensure equal

elapsed time for each pipetting step, without interruption. Duplication of

all standards and specimens, although not required, is recommended. To

avoid cross-contamination, change pipette tips between additions of each

standard level, between sample additions, and between reagent additions.

Also, use separate reservoirs for each reagent.

3. Incubation: To ensure accurate results, proper adhesion of plate

sealers during incubation steps is necessary. Do not allow wells to sit

uncovered for extended periods between incubation steps. Once

reagents have been added to the well strips, DO NOT let the strips DRY

at any time during the assay. Incubation time and temperature must be

observed.





Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com



- Washing: The wash procedure is critical. Complete removal of 4. liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Solution by aspirating or decanting and remove any drop of water and fingerprint on the bottom of the plate. Insufficient washing will result in poor precision and falsely elevated absorbance reading. When using an automated plate washer, adding a 30 second soak period following the addition of wash buffer, and/or rotating the plate 180 degrees between wash steps may improve assay precision.
- 5. Controlling of reaction time: Observe the change of color after adding TMB Substrate (e.g. observation once every 10 mins), TMB Substrate should change from colorless or light blue to gradations of blue. If the color is too deep, add Stop Solution in advance to avoid excessively strong reaction which will result in inaccurate absorbance reading.
- 6. TMB Substrate is easily contaminated. TMB Substrate should remain colorless or light blue until added to the plate. Please protect it from light.
- 7 Stop Solution should be added to the plate in the same order as the TMB Substrate. The color developed in the wells will turn from blue to yellow upon addition of the Stop Solution. Wells that are green in color indicate that Stop Solution has not mixed thoroughly with TMB Substrate.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Calculation of Results

BIOMATIK

Using the professional soft "Curve Expert" to make a standard

curve is recommended, which can be downloaded from our web.

Average the duplicate readings for each standard and sample and

subtract the average optical density of Blank.

Create a standard curve by reducing the data using computer software

capable of generating a four parameter logistic (4-PL) curve-fit. As an

alternative, construct a standard curve by plotting the mean absorbance

for each standard on the x-axis against the concentration on the y-axis

and draw a best fit curve through the points on the graph. The data may

be linearized by plotting the log of the Alb concentrations versus the log

of the O.D. and the best fit line can be determined by regression analysis.

This procedure will produce an adequate but less precise fit of the data.

If samples have been diluted, the concentration read from the standard

curve must be multiplied by the dilution factor.