



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com
https://www.biomatik.com

User Manual

Catalog Number: EKC34103

Product Name: Human Ig gamma-3 chain C region(IGHG3) ELISA kit

Detection Range: 2.7 ng/mL-2000 ng/mL

Intended Use: For quantitative determination of Human Ig gamma-3 chain C

region (IGHG3) concentrations in serum, plasma.

Precautions: For research use only. Not for use in diagnostic procedures.

Manual Version: 202301V1

Storage:

Unopened kit	6 months when stored at 2 - 8°C.
Opened Kit	May be stored up to 1 month at 2 - 8°C. Keep it in sealed aluminum foil bag and avoid moisture.

The product manual may be updated as a result of continuous improvements.

Always refer to the hard copy manual included in the kit for your experiment.





Email: info@biomatik.com https://www.biomatik.com

Kit Components

Reagents	Quantity		
Assay plate (12 x 8 coated Microwells)	1(96 wells)		
Standard (Freeze dried)	2		
HRP-conjugate (100 x concentrate)	1 x 60 µl		
HRP-conjugate Diluent	1 x 20 ml		
Sample Diluent	2 x 20 ml		
Wash Buffer (25 x concentrate)	1 x 20 ml		
TMB Substrate	1 x 10 ml		
Stop Solution	1 x 10 ml		
Adhesive Strip (For 96 wells)	4		
Instruction manual	1		

Working Principle

This assay employs the competitive inhibition enzyme immunoassay technique. The microtiter plate provided in this kit has been pre-coated with IGHG3. Standards or samples are added to the appropriate microtiter plate wells with Horseradish Peroxidase (HRP) conjugated antibody specific to IGHG3. The competitive inhibition reaction is launched between with pre-coated IGHG3 and IGHG3 in samples. A substrate solution is added to the wells and the color develops in opposite to the amount of IGHG3 in the sample. Measure the intensity of the color when the color stopped developing.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Sensitivity

BIOMATIK

The minimum detectable dose of human IGHG3 is less than 2.7 ng/mL.

The sensitivity of this assay, or Lower Limit of Detection (LLD) was

defined as the lowest human IGHG3 concentration that could be

differentiated from zero. It was determined by the mean O.D value of 20

replicates of the zero standard added by their three standard deviations.

Specificity

This assay has high sensitivity and excellent specificity. No significant

cross-reactivity or interference between Human IGHG3 and analogues

was observed. Limited by current skills and knowledge, it is impossible

for us to complete the cross-reactivity detection between Human IGHG3

and all the analogues, therefore, cross reaction may still exist.

Precision

Intra-assay Precision (Precision within an assay): CV%<8%

Three samples of known concentration were tested twenty times on one

plate to assess.

Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to

assess.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

BIOMATIK

Other Supplies Required

Microplate reader capable of measuring absorbance at 450 nm,

with the correction wavelength set at 540 nm or 570 nm.

An incubator which can provide stable incubation conditions up to

37°C±0.5°C.

Squirt bottle, manifold dispenser, or automated microplate washer.

Absorbent paper for blotting the microtiter plate.

100ml and 500ml graduated cylinders.

Deionized or distilled water.

Pipettes and pipette tips.

Test tubes for dilution.

Sample Collection & Storage

Serum Use a serum separator tube (SST) and allow samples to

clot for two hours at room temperature or overnight at 4°C before

centrifugation at 1000 ×g for 15 mins. Remove serum and assay

immediately or aliquot and store samples at -20°C or -80°C. Avoid

repeated freeze-thaw cycles.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Plasma Collect plasma using EDTA or heparin as an anticoagulant.

Within 30 mins after collecting samples, centrifuge samples at 1000 x g, 2

- 8°C, for 15 mins. Assay immediately or aliquot and store samples at -

20°C or -80°C. Avoid repeated freeze-thaw cycles.

BIOMATIK

Sample Preparation

Recommend to dilute the serum or plasma samples with Sample

Diluent(1:2000) before test. The suggested 2000-fold dilution can be

achieved by adding 2µl sample to 98µl of Sample Diluent. Complete the

2000-fold dilution by adding 6µl of this solution to 234µl of Sample Diluent.

The recommended dilution factor is for reference only. The optimal

dilution factor should be determined by users according to their particular

experiments.

Notes:

1. Biomatik is only responsible for the kit itself, not for the samples

consumed during the assay. The user need to calculate the possible

amount of the samples to be used in the whole test. Please reserve

sufficient samples in advance.

Email: info@biomatik.com https://www.biomatik.com

- 2. Samples to be used within 5 days may be stored at 2-8°C, otherwise, samples must be stored at -20°C (≤1month) or -80°C (≤2month) to avoid contamination and loss of bioactivity.
- 3. Grossly hemolyzed samples are not suitable to use.
- 4. It would be necessary to run a preliminary experiment for validation, if the samples are not indicated in the manual.
- 5. Please predict the concentration before assaying. If results were not within the range of the standard curve, users would need to estimate the optimal sample dilutions for their particular experiments.
- Tissue or cell extraction samples prepared by chemical lysis buffer may cause unexpected ELISA results due to the impacts of certain chemicals.
- 7. Considering the possibility of mismatch between antigen from other resources and antibody in our kits (e.g., antibody targets conformational epitope rather than linear epitope), some native or recombinant proteins from other manufacturers may not be recognized by our products.
- Due to factors including cell viability, cell number and sampling time, samples from cell culture supernatant may not be detected by the kit.
- Recommend to use fresh samples for the test. If you store samples for long time, protein degradation and denaturalization may occur in those samples and finally lead to wrong results.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Reagents Preparation

BIOMATIK

Kindly use graduated containers to prepare the reagent.

Please don't prepare the reagent directly in the Diluent vials in the

kit.

Bring all reagents to room temperature (18-25°C) before use for 30

mins.

Prepare fresh standard for each assay. Use within 4 hours and

discard after use

Making serial dilution in the wells directly is not permitted.

Distilled water is recommended to be used to make the preparation

for reagents. Contaminated water or container for reagent preparation will

influence the detection result.

HRP-conjugate (1x) - Centrifuge the vial before opening.

HRP-conjugate requires a 100-fold dilution. A suggested 100-fold

dilution is 10 µl HRP-conjugate + 990 µl HRP-conjugate Diluent.

Wash Buffer(1x)- If crystals have formed in the concentrate, warm up to

room temperature and mix gently until the crystals have completely

dissolved. Dilute 20 ml of Wash Buffer Concentrate (25 x) into deionized

or distilled water to prepare 500 ml of Wash Buffer (1 x).



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515


Email: info@biomatik.com https://www.biomatik.com

Standard

Centrifuge the standard vial at 6000-10000rpm for 30s.

Reconstitute the Standard with 1.0 ml Sample Diluent. Do not substitute other diluents. This reconstitution produces a stock solution of 2000 ng/mL. Mix the standard to ensure complete reconstitution and allow the standard to sit for a minimum of 15 mins with gentle agitation prior to making dilutions.

Pipette 200 µl Sample Diluent into each tube (S0-S6). Use the stock solution to produce a 3-fold dilution series (below). Mix each tube thoroughly before the next transfer. The undiluted Standard serves as the high standard (2000 ng/mL). Sample Diluent serves as the zero standard (0 ng/mL).

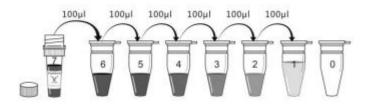

Tube	S7	S6	S5	\$4	S3	S2	S1	S0
ng/ml	2000	667	222	74	25	8.2	2.7	0

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Key Notes

Do not mix or substitute reagents with those from other lots or

sources

If samples generate higher values than the highest standard, dilute

the samples with Sample Diluent and repeat the assay.

Any variation in Sample Diluent, operator, pipetting technique,

washing technique, incubation time or temperature, and kit age can

cause variation in binding.

BIOMATIK

This assay is designed to eliminate interference by soluble

receptors, binding proteins, and other factors present in biological

samples. Until all factors have been tested in the Immunoassay, the

possibility of interference cannot be excluded.

Precautions

The Stop Solution provided with this kit is an acid solution. Wear eye, hand, face,

and clothing protection when using this material.

Assav Procedures

Bring all reagents and samples to room temperature before use.

Centrifuge the sample again after thawing before the assay. It is

recommended to assay all samples and standards in duplicate.



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

- Prepare all reagents, working standards, and samples as directed in the previous sections.
- Determine the number of wells to be used and put any remaining wells and the desiccant back into the pouch and seal the ziploc, store unused wells at 4°C.
- 3. Set a Blank well without any solution.
- 4. Add 50µl **standard** or **sample** per well. Add 50µl **HRP-conjugate** (1x) to each well immediately (not to Blank well). Mix well with the pipette or shake the plate gently for 60 seconds. A plate layout is provided to record standards and samples assayed.
- 5. Seal plate with the adhesive strip. Incubate at 37°C for 40 mins .
- 6. Aspirate each well and wash, repeating the process four times for a total of five washes. Wash by filling each well with **Wash Buffer** (200μl) using a squirt bottle, multi-channel pipette, manifold dispenser, or autowasher, and let it stand for 2 mins, complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.
- 7. Add $90\mu l$ TMB Substrate to each well. Incubate at $37^{\circ}C$ for 20 mins. Protect from light.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

8. Add 50µl Stop Solution to each well, gently tap the plate to ensure

9. Determine the optical density of each well within 5 mins, using a

microplate reader set to 450 nm. If wavelength correction is available, set

to 540 nm or 570 nm. Subtract readings at 540 nm or 570 nm from the

readings at 450 nm. This subtraction will correct for optical imperfections

in the plate. Readings made directly at 450 nm without correction may be

higher and less accurate.

thorough mixing.

BIOMATIK

*Samples may require dilution. Please refer to Sample Preparation

section.

Notes:

1. The final experimental results will be closely related to validity of

the products, operation skills of the end users and the experimental

environments.

2 Samples or reagents addition: Please use the freshly prepared

Standard. Please carefully add samples to wells and mix gently to avoid

foaming. Do not touch the well wall as possible. For each step in the

procedure, total dispensing time for addition of reagents or samples to

the assay plate should not exceed 10 mins. This will ensure equal





Email: info@biomatik.com

BIOMATIK

elapsed time for each pipetting step, without interruption. Duplication of all standards and specimens, although not required, is recommended. To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.

- 3. Incubation: To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary. Do not allow wells to sit uncovered for extended periods between incubation steps. Once reagents have been added to the well strips, DO NOT let the strips DRY at any time during the assay. Incubation time and temperature must be observed.
- 4. Washing: The wash procedure is critical. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Solution by aspirating or decanting and remove any drop of water and fingerprint on the bottom of the plate. Insufficient washing will result in poor precision and falsely elevated absorbance reading. When using an automated plate washer, adding a 30 second soak period following the addition of wash buffer, and/or rotating the plate 180 degrees between wash steps may improve assay precision.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

5. Controlling of reaction time: Observe the change of color after

adding TMB Substrate (e.g. observation once every 10 mins), TMB

Substrate should change from colorless or light blue to gradations of blue.

If the color is too deep, add Stop Solution in advance to avoid excessively

strong reaction which will result in inaccurate absorbance reading.

TMB Substrate is easily contaminated. TMB Substrate should 6.

remain colorless or light blue until added to the plate. Please protect it

from light.

7. Stop Solution should be added to the plate in the same order as

the TMB Substrate. The color developed in the wells will turn from blue to

yellow upon addition of the Stop Solution. Wells that are green in color

indicate that Stop Solution has not mixed thoroughly with TMB Substrate.

Calculation of Results

BIOMATIK

Using the professional soft "Curve Expert" to make a standard

curve is recommended, which can be downloaded from our web.

Average the duplicate readings for each standard and sample and

subtract the average optical density of Blank.





Email: info@biomatik.com

https://www.biomatik.com

Create a standard curve by reducing the data using computer software

capable of generating a four parameter logistic (4-PL) curve-fit. As an

alternative, construct a standard curve by plotting the mean absorbance

for each standard on the x-axis against the concentration on the y-axis

and draw a best fit curve through the points on the graph. The data may

be linearized by plotting the log of the IGHG3 concentrations versus the

log of the O.D. and the best fit line can be determined by regression

analysis. This procedure will produce an adequate but less precise fit of

the data.

If samples have been diluted, the concentration read from the standard

curve must be multiplied by the dilution factor.