

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

User Manual

Catalog Number: EKC34395

Product Name: Human leucine-rich alpha-2 glycoprotein 1(LRG1) ELISA

Kit

Detection Range: 0.156 µg/ml-40 µg/ml

Intended Use: For quantitative determination of human leucine-rich alpha-2

glycoprotein 1 (LRG1) concentrations in serum, urine, tissue homogenates.

Precautions: For research use only. Not for use in diagnostic procedures.

Manual Version: 202301V1

Storage:

Unopened kit	6 months when stored at 2 - 8°C.
Opened Kit	May be stored up to 1 month at 2 - 8°C. Keep it in sealed
Opened Kit	aluminum foil bag and avoid moisture.

The product manual may be updated as a result of continuous improvements.

Always refer to the hard copy manual included in the kit for your experiment.



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

Kit Components

Reagents	Quantity
Assay Plate (12 x 8 coated Microwells)	1(96 wells)
Standard	2 x 250 µl
Antibody (100 x concentrate)	1 x 60 µl
HRP-conjugate (100 x concentrate)	1 x 120 µl
Antibody Diluent	1 x 10 ml
HRP-conjugate Diluent	1 x 20 ml
Sample Diluent	2 x 20 ml
Wash Buffer (25 x concentrate)	1 x 20 ml
TMB Substrate	1 x 10 ml
Stop Solution	1 x 10 ml
Adhesive Strip (For 96 wells)	4
Instruction Manual	1

Working Principle

This assay employs the competitive inhibition enzyme immunoassay technique. The microtiter plate provided in this kit has been pre-coated with an antigen. Standards or samples are added to the appropriate microtiter plate wells with antibody specific for LRG1 and Horseradish Peroxidase (HRP) conjugated goat-anti-rabbit antibody. The competitive inhibition reaction is launched between with pre-coated LRG1 and LRG1 in samples. A substrate solution is added to the wells and the color develops in opposite to the amount of LRG1 in the sample. At last, measure the intensity of color after stopping color development.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Sensitivity

The minimum detectable dose of human LRG1 is less than 0.156 µg/ml.

The sensitivity of this assay, or Lower Limit of Detection (LLD) was

defined as the lowest human LRG1 concentration that could be

differentiated from zero.

BIOMATIK

Specificity

This assay has high sensitivity and excellent specificity. No significant

cross-reactivity or interference between human LRG1 and analogues was

observed. Limited by current skills and knowledge, it is impossible for us

to complete the cross-reactivity detection between human LRG1 and all

the analogues, therefore, cross reaction may still exist.

Precision

Intra-assay Precision (Precision within an assay): CV%<12%

Three samples of known concentration were tested twenty times on one

plate to assess.

Inter-assay Precision (Precision between assays): CV%<15%

Three samples of known concentration were tested in twenty assays to

assess.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

BIOMATIK

Other Supplies Required

Microplate reader capable of measuring absorbance at 450 nm,

with the correction wavelength set at 600 nm or 630 nm.

An incubator which can provide stable incubation conditions up to

37°C±0.5°C.

Squirt bottle, manifold dispenser, or automated microplate washer.

Absorbent paper for blotting the microtiter plate.

100ml and 500ml graduated cylinders.

Deionized or distilled water.

Pipettes and pipette tips.

Test tubes for dilution.

Sample Collection & Storage

Serum Use a serum separator tube (SST) and allow samples to

clot for two hours at room temperature or overnight at 4°C before

centrifugation at 1000 ×g for 15 mins. Remove serum and assay

immediately or aliquot and store samples at -20°C or -80°C. Avoid

repeated freeze-thaw cycles.



- Tissue Homogenates 100mg tissue was rinsed with 1X PBS, homogenized in 1 ml of 1X PBS and stored overnight at -20°C. After two freeze-thaw cycles were performed to break the cell membranes, the homogenates were centrifuged for 5 minutes at 5000 x g, 2 - 8°C. The supernate was removed and assayed immediately. Alternatively, aliquot and store samples at -20°C or -80°C. Centrifuge the sample again after thawing before the assay. Avoid repeated freeze-thaw cycles.
- **Urine** Use a sterile container to collect urine samples. Remove any particulates by centrifugation for 15 minutes at 1000x g, 2 - 8°C and assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles. Centrifuge again before assaying to remove any additional precipitates that may appear after storage.

Notes:

- 1. Biomatik is only responsible for the kit itself, not for the samples consumed during the assay. The user need to calculate the possible amount of the samples to be used in the whole test. Please reserve sufficient samples in advance.
- 2. Samples to be used within 5 days may be stored at 2-8°C, otherwise. samples must be stored at -20°C (≤1month) or -80°C (≤2month) to avoid contamination and loss of bioactivity.

BIOMATIK

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

3. Grossly hemolyzed samples are not suitable to use.

4. It would be necessary to run a preliminary experiment for validation, if

the samples are not indicated in the manual.

5. Please predict the concentration before assaying. If results were not

within the range of the standard curve, users would need to estimate the

optimal sample dilutions for their particular experiments.

6. Tissue or cell extraction samples prepared by chemical lysis buffer

may cause unexpected ELISA results due to the impacts of certain

chemicals.

7. Considering the possibility of mismatch between antigen from other

resources and antibody in our kits (e.g., antibody targets conformational

epitope rather than linear epitope), some native or recombinant proteins

from other manufacturers may not be recognized by our products.

8. Due to factors including cell viability, cell number and sampling time,

samples from cell culture supernatant may not be detected by the kit.

9. Recommend to use fresh samples for the test. If you store samples for

long time, protein degradation and denaturalization may occur in those

samples and finally lead to wrong results.

Reagents Preparation

- Kindly use graduated containers to prepare the reagent. Please don't prepare the reagent directly in the Diluent vials provided in the kit.
- Bring all reagents to room temperature (18-25°C) before use for 30 mins.
- Prepare fresh standard for each assay. Use within 4 hours and discard after use.
- Making serial dilution in the wells directly is not permitted.
- Distilled Water is recommended. Contaminated water or container for reagents preparation will affect the test result.

Antibody (1x) - Centrifuge the vial before opening.

Antibody requires a 100-fold dilution. A suggested 100-fold dilution is 10 µl of Antibody + 990 µl of Antibody Diluent.

HRP-conjugate (1x) - Centrifuge the vial before opening.

HRP-conjugate requires a 100-fold dilution. A suggested 100-fold dilution is 10 µl of HRP-conjugate + 990 µl of HRP-conjugate Diluent.





Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515


Email: info@biomatik.com https://www.biomatik.com

Wash Buffer (1x)- If crystals have formed in the concentrate, warm them up to room temperature and mix them gently until they get completely dissolved. Dilute 20 ml of Wash Buffer Concentrate (25 x) into deionized or distilled water to prepare 500 ml of Wash Buffer (1 x).

Standard

Centrifuge the standard vial at 6000-10000rpm for 30s before opening. Allow the stock standard solution of 40 µg/ml to sit for a minimum of 30 mins to warm to room temperature prior to making dilutions.

Pipette 150 µl Sample Diluent into each tube (S0-S4). Use the stock solution to produce a 4-fold dilution series (below). Mix each tube thoroughly before the next transfer. The undiluted Standard serves as the high standard (40 μg/ml). **Sample Diluent** serves as the zero standard (0 μg/ml).

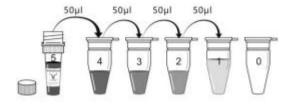





Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140. (877) 221-3515

Email: info@biomatik.com

Tube	S5	S4	S3	S2	S1	S0
µg/ml	40	10	2.5	0.625	0.156	0

Key Notes

- Do not mix or substitute reagents with those from other lots or sources
- If samples generate higher values than the highest standard, dilute the samples with Sample Diluent and repeat the assay.
- Any variation in Sample Diluent, operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding.
- This assay is designed to eliminate interference by soluble receptors, binding proteins, and other factors present in biological samples. Until all factors have been tested in the Immunoassay, the possibility of interference cannot be excluded.

Precautions

The Stop Solution provided with this kit is an acid solution. Wear eye, hand, face, and clothing protection when using this material.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Assay Procedures

BIOMATIK

Bring all reagents and samples to room temperature before use.

Centrifuge the sample again after thawing before the assay. It is

recommended to assay all samples and standards in duplicate.

1. Prepare all reagents, working standards, and samples as directed in

the previous sections.

2. Determine the number of wells to be used and put remaining wells and the

desiccant back into the pouch and seal the ziploc, store unused wells at 4°C.

3. Set a **Blank** well without any solution.

4. Add 50µl Standard or Sample per well. Add 50µl Antibody (1x) to each

well immediately (not to Blank well). Mix well with the pipette or shake the

plate gently for 60 seconds. A plate layout is provided to record standards

and samples assayed.

5. Seal plate with adhesive strip. Incubate at 37°C for 30 mins.

6. Aspirate each well and wash, repeating the process two times for a total of

three washes. Wash by filling each well with Wash Buffer (250µl) using a

squirt bottle, multi-channel pipette, manifold dispenser, or autowasher, and

let it stand for 2 mins, complete removal of liquid at each step is essential to

good performance. After the last wash, remove any remaining Wash Buffer

by aspirating or decanting. Invert the plate and blot it against clean paper

towels

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

BIOMATIK

7. Add 100µl HRP-conjugate (1x) to each well immediately (not to Blank well). Seal plate with adhesive strip. Incubate at 37°C for 30 mins.

8. Repeat the aspiration/wash process for five times as in step 6.

9. Add 90µl TMB Substrate to each well. Incubate at 37°C for 20 mins.

Protect from light.

10. Add 50µl Stop Solution to each well, gently tap the plate to ensure

thorough mixing.

11. Determine the optical density of each well within 5 mins, using a

microplate reader set to 450 nm. If wavelength correction is available, set to

540 nm or 570 nm. Subtract readings at 540 nm or 570 nm from the readings

at 450 nm. This subtraction will correct for optical imperfections in the plate.

Readings made directly at 450 nm without correction may be higher and less

accurate

Notes:

1. The final experimental results will be closely related to validity of

products, operation skills of end users and the experimental

environments

2 Samples or reagents addition: Please use the freshly prepared

Standard. Please carefully add samples to wells and mix gently to avoid

foaming. Do not touch the well wall as possible. For each step in the



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

BIOMATIK

procedure, total dispensing time for addition of reagents or samples to the assay plate should not exceed 10 mins. This will ensure equal elapsed time for each pipetting step, without interruption. Duplication of all standards and specimens, although not required, is recommended. To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.

- 3 Incubation: To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary. Do not allow wells to sit uncovered for extended periods between incubation steps. Once reagents have been added to the well strips, DO NOT let the strips DRY at any time during the assay. Incubation time and temperature must be observed.
- Washing: The wash procedure is critical. Complete removal of 4. liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Solution by aspirating or decanting and remove any drop of water and fingerprint on the bottom of the plate. Insufficient washing will result in poor precision and falsely elevated absorbance reading. When using an automated plate washer, adding a 30 second soak period following the addition of wash buffer, and/or rotating the plate 180 degrees between wash steps may improve assay

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

precision.

BIOMATIK

Controlling of reaction time: Observe the change of color after 5

adding TMB Substrate (e.g. observation once every 10 mins), TMB

Substrate should change from colorless or light blue to gradations of blue.

If the color is too deep, add Stop Solution in advance to avoid excessively

strong reaction which will result in inaccurate absorbance reading.

TMB Substrate is easily contaminated. TMB Substrate should 6.

remain colorless or light blue until added to the plate. Please protect it

from light.

7. Stop Solution should be added to the plate in the same order as

the TMB Substrate. The color developed in the wells will turn from blue to

yellow upon addition of the Stop Solution. Wells that are green in color

indicate that Stop Solution has not mixed thoroughly with TMB Substrate.

Calculation of Results

Using the professional soft "Curve Expert" to make a standard

curve is recommended, which can be downloaded from our web.

Average the duplicate readings for each standard and sample and

subtract the average optical density of Blank.





Tel: (519) 489-7195, (800) 836-8089

Fax: (519) 231-0140, (877) 221-3515 Email: info@biomatik.com

https://www.biomatik.com

Create a standard curve by reducing the data using computer software

capable of generating a four parameter logistic (4-PL) curve-fit. As an

alternative, construct a standard curve by plotting the mean absorbance

for each standard on the x-axis against the concentration on the y-axis

and draw a best fit curve through the points on the graph. The data may

be linearized by plotting the log of the LRG1 concentrations versus the

log of the O.D. and the best fit line can be determined by regression

analysis. This procedure will produce an adequate but less precise fit of

the data.

If samples have been diluted, the concentration read from the standard

curve must be multiplied by the dilution factor.