

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com
https://www.biomatik.com

User Manual

Catalog Number: EKC34890

Product Name: Human orosomucoid 2 (ORM2) ELISA Kit

Detection Range: 78 ng/ml-5000 ng/ml

Intended Use: For quantitative determination of human orosomucoid 2 (ORM2)

concentrations in serum, plasma and tissue homogenates.

Precautions: For research use only. Not for use in diagnostic procedures.

Manual Version: 202301V1

Storage:

Unopened kit	6 months when stored at 2 - 8°C.
Opened Kit	May be stored up to 1 month at 2 - 8°C. Keep it in sealed aluminum foil bag and avoid moisture.

The product manual may be updated as a result of continuous improvements.

Always refer to the hard copy manual included in the kit for your experiment.



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

Kit Components

Reagents	Quantity		
Assay Plate (12 x 8 coated Microwells)	1(96 wells)		
Standard (Freeze dried)	2		
HRP-conjugate (100 x concentrate)	1 x 60 µl		
HRP-conjugate Diluent	1 x 10 ml		
Sample Diluent	2 x 20 ml		
Wash Buffer (25 x concentrate)	1 x 20 ml		
TMB Substrate	1 x 10 ml		
Stop Solution	1 x 10 ml		
Adhesive Strip (For 96 wells)	4		
Instruction Manual	1		

Working Principle

This assay employs the competitive inhibition enzyme immunoassay technique. The microtiter plate provided in this kit has been pre-coated with ORM2. Standards or samples are added to the appropriate microtiter plate wells with Horseradish Peroxidase (HRP) conjugated antibody preparation specific for ORM2. The competitive inhibition reaction is launched between with pre-coated ORM2 and ORM2 in samples. A substrate solution is added to the wells and the color develops in opposite to the amount of ORM2 in the samples. At last, measure the intensity of color after stopping color development.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140. (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Sensitivity

BIOMATIK

The minimum detectable dose of human ORM2 is less than 19.5 ng/ml.

The sensitivity of this assay, or Lower Limit of Detection (LLD) was

defined as the lowest human ORM2 concentration that could be

differentiated from zero. It was determined the mean O.D value of 20

replicates of the zero standard added by their three standard deviations.

Specificity

This assay has high sensitivity and excellent specificity. No significant

cross-reactivity or interference between human ORM2 and analogues

was observed. Limited by current skills and knowledge, it is impossible

for us to complete the cross-reactivity detection between human ORM2

and all the analogues, therefore, cross reaction may still exist.

Precision

Intra-assay Precision (Precision within an assay): CV%<8%

Three samples of known concentration were tested twenty times on one

plate to assess.





Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com
https://www.biomatik.com

Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to assess.

Other Supplies Required

- Microplate reader capable of measuring absorbance at 450 nm,
 with the correction wavelength set at 600 nm or 630 nm.
- An incubator which can provide stable incubation conditions up to 37°C±0.5°C.
- Squirt bottle, manifold dispenser, or automated microplate washer.
- Absorbent paper for blotting the microtiter plate.
- 100ml and 500ml graduated cylinders.
- Deionized or distilled water.
- Pipettes and pipette tips.
- Test tubes for dilution.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140. (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Sample Collection & Storage

(2) BIOMATIK

Serum Use a serum separator tube (SST) and allow samples to

clot for two hours at room temperature or overnight at 4°C before

centrifugation for 15 minutes at 1000 xg. Remove serum and assay

immediately or aliquot and store samples at -20°C or -80°C. Avoid

repeated freeze-thaw cycles.

Plasma Collect plasma using EDTA, or heparin as an

anticoagulant. Centrifuge for 15 minutes at 1000 ×g at 2-8°C within 30

minutes of collection. Assay immediately or aliquot and store samples at -

20°C or -80°C. Avoid repeated freeze-thaw cycles.

Tissue Homogenates 100mg tissue was rinsed with 1X PBS,

homogenized in 1 ml of 1X PBS and stored overnight at -20°C. After two

freeze-thaw cycles were performed to break the cell membranes, the

homogenates were centrifuged for 5 minutes at 5000 x g, 2 - 8°C. The

supernate was removed and assayed immediately. Alternatively, aliquot

and store samples at -20°C or -80°C. Centrifuge the sample again after

thawing before the assay. Avoid repeated freeze-thaw cycles.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140. (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Sample Preparation

BIOMATIK

Recommend to dilute the serum or plasma samples with Sample

Diluent(1:500) before test. The suggested 500-fold dilution can be

achieved by adding 5µl sample to 95µl of Sample Diluent first, then

complete the 500-fold dilution by adding 10µl of this solution to 240µl of

Sample Diluent. The recommended dilution factor is for reference only.

The optimal dilution factor should be determined by users according to

their particular experiments.

Notes:

1. Biomatik is only responsible for the kit itself, not for the samples

consumed during the assay. The user need to calculate the possible

amount of the samples to be used in the whole test. Please reserve

sufficient samples in advance.

2. Samples to be used within 5 days may be stored at 2-8°C, otherwise,

samples must be stored at -20°C (≤1month) or -80°C (≤2month) to avoid

contamination and loss of bioactivity.

3. Grossly hemolyzed samples are not suitable to use.

4. It would be necessary to run a preliminary experiment for validation, if

the samples are not indicated in the manual.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

BIOMATIK

5. Please predict the concentration before assaying. If results were not within the range of the standard curve, users would need to estimate the optimal sample dilutions for their particular experiments.

6. Tissue or cell extraction samples prepared by chemical lysis buffer may cause unexpected ELISA results due to the impacts of certain chemicals.

7. Considering the possibility of mismatch between antigen from other resources and antibody in our kits (e.g., antibody targets conformational epitope rather than linear epitope), some native or recombinant proteins from other manufacturers may not be recognized by our products.

8. Due to factors including cell viability, cell number and sampling time, samples from cell culture supernatant may not be detected by the kit.

9. Recommend to use fresh samples for the test. If you store samples for long time, protein degradation and denaturalization may occur in those samples and finally lead to wrong results.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Reagents Preparation

BIOMATIK

Kindly use graduated containers to prepare the reagent.

Please don't prepare the reagent directly in the Diluent vials

provided in the kit.

Bring all reagents to room temperature (18-25°C) before use for 30

mins.

Prepare fresh standard for each assay. Use within 4 hours and

discard after use

Making serial dilution in the wells directly is not permitted.

Please carefully reconstitute Standards according to the instruction,

and avoid foaming and mix gently until the crystals have completely

dissolved. To minimize imprecision caused by pipetting, use small

volumes and ensure that pipettors are calibrated. It is recommended to

suck more than 10µl for once pipetting.

Distilled water is recommended to be used to make the preparation

for reagents. Contaminated water or container for reagent preparation will

influence the detection result

HRP-conjugate (1x) - Centrifuge the vial before opening.

HRP-conjugate requires a 100-fold dilution. A suggested 100-fold dilution is

10 ul of HRP-conjugate + 990 ul of HRP-conjugate Diluent.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Wash Buffer (1x)- If crystals have formed in the concentrate, warm them up

to room temperature and mix them gently until they get completely dissolved.

Dilute 20 ml of Wash Buffer Concentrate (25 x) into deionized or distilled

water to prepare 500 ml of Wash Buffer (1 x).

BIOMATIK

Standard

Centrifuge the standard vial at 6000-10000rpm for 30s.

Reconstitute the Standard with 1.0 ml of Sample Diluent. Do not substitute

other diluents. This reconstitution produces a stock solution of 5000 ng/ml.

Mix the standard to ensure complete reconstitution and allow the standard to

sit for a minimum of 15 mins with gentle agitation prior to making dilutions.

Pipette 150 µl Sample Diluent into each tube (S0-S6). Use the stock

solution to produce a 2-fold dilution series (below). Mix each tube thoroughly


before the next transfer. The undiluted Standard serves as the high standard

(5000 ng/ml). Sample Diluent serves as the zero standard (0 ng/ml).



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140. (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

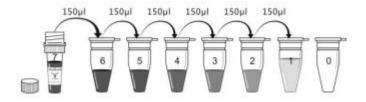

Tube	S7	S6	S5	S4	S3	S2	S1	SO
ng/ml	5000	2500	1250	625	312	156	78	0

Key Notes

- Do not mix or substitute reagents with those from other lots or sources.
- If samples generate higher values than the highest standard, dilute the samples with Sample Diluent and repeat the assay.
- Any variation in Sample Diluent, operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding.
- This assay is designed to eliminate interference by soluble receptors, binding proteins, and other factors present in biological samples. Until all factors have been tested in the Immunoassay, the possibility of interference cannot be excluded.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Precautions

BIOMATIK

The Stop Solution provided with this kit is an acid solution. Wear eve. hand, face.

and clothing protection when using this material.

Assay Procedures

Bring all reagents and samples to room temperature before use.

Centrifuge the sample again after thawing before the assay. It is

recommended to assay all samples and standards in duplicate.

1. Prepare all reagents, working standards, and samples as directed in

the previous sections.

2. Determine the number of wells to be used and put remaining wells and the

desiccant back into the pouch and seal the ziploc, store unused wells at 4°C.

3. Set a **Blank** well without any solution.

4. Add 50µl Standard or Sample per well.

5. Add 50µl HRP-conjugate(1x) to each well at once (Not to the Blank well).

Seal plate with adhesive strip. Mix well and then incubate at 37°C for 60 mins.

A plate layout is provided to record standards and samples assayed.

6. Aspirate each well and wash, repeat the process four times for a total of

five washes. Wash by filling each well with Wash Buffer (200µl) using a

squirt bottle, multi-channel pipette, manifold dispenser, or autowasher, and

let it stand for 2 mins. Complete removal of liquid at each step is essential to

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

good performance. After the last wash, remove any remaining Wash Buffer

by aspirating or decanting. Invert the plate and blot it against clean paper

towels

7. Add 90µl **TMBSubstrate** to each well and mix well. Cover the microtiter

plate adhesive strip. Incubate at 37°C for 20 mins. Protect from light.

8. Add 50µl Stop Solution to each well, gently tap the plate to ensure

thorough mixing.

BIOMATIK

9. Determine the optical density of each well within 5 mins, using a

microplate reader set to 450 nm. If wavelength correction is available, set to

540 nm or 570 nm. Subtract readings at 540 nm or 570 nm from the readings

at 450 nm. This subtraction will correct for optical imperfections in the plate.

Readings made directly at 450 nm without correction may be higher and less

accurate.

*Samples may require dilution. Please refer to Sample Preparation.

Notes:

1. The final experimental results will be closely related to validity of

products, operation skills of end users and the experimental

environments.





Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com
https://www.biomatik.com

2. Samples or reagents addition: Please use the freshly prepared Standard. Please carefully add samples to wells and mix gently to avoid foaming. Do not touch the well wall as possible. For each step in the procedure, total dispensing time for addition of reagents or samples to the assay plate should not exceed 10 mins. This will ensure equal elapsed time for each pipetting step, without interruption. Duplication of all standards and specimens, although not required, is recommended. To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.

- 3. Incubation: To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary. Do not allow wells to sit uncovered for extended periods between incubation steps. Once you add reagents to the well strips, DO NOT let the strips DRY at any time during the assay. Incubation time and temperature must be observed.
- 4. Washing: The wash procedure is critical. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Solution by aspirating or decanting and remove any drop of water and fingerprint on the bottom of the plate. Insufficient washing will result in poor precision and falsely elevated absorbance reading. When using an automated plate washer, adding a





Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140. (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

30 second soak period following the addition of wash buffer, and/or rotating the plate 180 degrees between wash steps may improve assay precision.

- 5. Controlling of reaction time: Observe the change of color after adding TMB Substrate (e.g. observation once every 10 mins), TMB Substrate should change from colorless or light blue to gradations of blue. If the color is too deep, add Stop Solution in advance to avoid excessively strong reaction which will result in inaccurate absorbance reading.
- TMB Substrate is easily contaminated. TMB Substrate should remain colorless or light blue until added to the plate. Please protect it from light.
- 7. Stop Solution should be added to the plate in the same order as the TMB Substrate. The color developed in the wells will turn from blue to yellow upon addition of the Stop Solution. Wells that are green in color indicate that Stop Solution has not mixed thoroughly with TMB Substrate.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Calculation of Results

BIOMATIK

Using the professional soft "Curve Expert" to make a standard

curve is recommended, which can be downloaded from our web.

Average the duplicate readings for each standard and sample and

subtract the average optical density of Blank.

Create a standard curve by reducing the data using computer software

capable of generating a four parameter logistic (4-PL) curve-fit. As an

alternative, construct a standard curve by plotting the mean absorbance

for each standard on the x-axis against the concentration on the y-axis

and draw a best fit curve through the points on the graph. The data may

be linearized by plotting the log of the ORM2 concentrations versus the

log of the O.D. and the best fit line can be determined by regression

analysis. This procedure will produce an adequate but less precise fit of

the data.

If samples have been diluted, the concentration read from the standard

curve must be multiplied by the dilution factor.