

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

# **User Manual**

Catalog Number: EKC39123

Product Name: Rat glycated hemoglobin A1c (GHbA1c) ELISA Kit

Detection Range: 0.019 ng/ml-80 ng/ml

Intended Use: For quantitative determination of rat glycated hemoglobin A1c

(GHbA1c) concentrations in lysate for RBC.

Precautions: For research use only. Not for use in diagnostic procedures.

Manual Version: 202301V1

### Storage:

Unopened kit	6 months when stored at 2 - 8°C.
Opened Kit	May be stored up to 1 month at 2 - 8°C. Keep it in sealed aluminum foil bag and avoid moisture.

The product manual may be updated as a result of continuous improvements.

Always refer to the hard copy manual included in the kit for your experiment.



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

### Kit Components

Reagents	Quantity		
Assay Plate (12 x 8 coated Microwells)	1(96 wells)		
Standard (10 x concentrate)	1 x 200 µl		
Biotin-conjugate (100 x concentrate)	1 x 60 µl		
HRP-avdin (100 x concentrate)	1 x 120 µl		
Biotin-conjugate Diluent	1 x 10 ml		
HRP-avdin Diluent	1 x 20 ml		
Sample Diluent	2 x 20 ml		
Wash Buffer (25 x concentrate)	1 x 20 ml		
TMB Substrate	1 x 10 ml		
Stop Solution	1 x 10 ml		
Adhesive Strip (For 96 wells)	4		
Instruction Manual	1		

### **Working Principle**

This assay employs the direct competitive inhibition enzyme immunoassay technique. Antibody specific for GHbA1c has been precoated onto a microplate. Standards and samples are pipetted into the wells with biotin-conjugated GHbA1c. A competitive inhibition reaction is launched between GHbA1c (Standards or samples) and biotin-conjugated GHbA1c with the pre-coated GHbA1c antibody. After washing, avidin conjugated Horseradish Peroxidase (HRP) is added to the wells. Following a wash to remove any unbound reagent, a substrate solution is added to the wells and color develops in opposite to the amount of GHbA1c bound in the initial step. At last, measure the intensity of color after stopping color development.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Sensitivity

The minimum detectable dose of rat GHbA1c is less than 0.019 pg/ml.

The sensitivity of this assay, or Lower Limit of Detection (LLD) was

defined as the lowest rat GHbA1c concentration that could be

differentiated from zero.

**BIOMATIK** 

Specificity

This assay has high sensitivity and excellent specificity. No significant

cross-reactivity or interference between rat GHbA1c and analogues was

observed. Limited by current skills and knowledge, it is impossible for us

to complete the cross-reactivity detection between rat GHbA1c and all the

analogues, therefore, cross reaction may still exist.

Precision

Intra-assay Precision (Precision within an assay): CV%<8%

Three samples of known concentration were tested twenty times on one

plate to assess.

Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to

assess.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Other Supplies Required

**BIOMATIK** 

Microplate reader capable of measuring absorbance at 450 nm,

with the correction wavelength set at 600 nm or 630 nm.

An incubator which can provide stable incubation conditions up to

37°C±0.5°C.

Squirt bottle, manifold dispenser, or automated microplate washer.

Absorbent paper for blotting the microtiter plate.

100ml and 500ml graduated cylinders.

Deionized or distilled water.

Pipettes and pipette tips.

Test tubes for dilution.

Sample Collection & Storage

Lysate for RBC Collect a certain volume of anticoagulant blood

sample (e.g. 1ml). Centrifuge at 10000rpm for 10min at 4°C to collect the

RBC. Remove supernate and add cold pure water to redissolve RBC to

the previous volume (e.g. 1ml). Then store at -20°C for about 30-60

minutes. After three freeze-thaw cycles to break up the cell membranes,

the cell lysates were centrifuged for 10 minutes at 10000rpm, 4°C, Collect

the supernatant. RBC lysates should be assayed immediately or

aliquotted and stored at -20°C. Centrifuge the sample again after thawing

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

before the assay. Avoid repeated freeze-thaw cycles.

Sample Preparation

**BIOMATIK** 

Recommend to dilute the lysate for RBC samples with Sample

Diluent(1:5000) before test. The suggested 5000-fold dilution can be

achieved by adding 2µl sample to 98µl of Sample Diluent. Complete the

5000-fold dilution by adding 3µl of this solution to 297µl of Sample Diluent.

The recommended dilution factor is for reference only. The optimal

dilution factor should be determined by users according to their particular

experiments.

Notes:

1. Biomatik is only responsible for the kit itself, not for the samples

consumed during the assay. The user need to calculate the possible

amount of the samples to be used in the whole test. Please reserve

sufficient samples in advance.

2. Samples to be used within 5 days may be stored at 2-8°C, otherwise,

samples must be stored at -20°C (≤1month) or -80°C (≤2month) to avoid

contamination and loss of bioactivity.

3. Grossly hemolyzed samples are not suitable to use.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

**BIOMATIK** 

4. It would be necessary to run a preliminary experiment for validation, if the samples are not indicated in the manual.

5. Please predict the concentration before assaying. If results were not within the range of the standard curve, users would need to estimate the optimal sample dilutions for their particular experiments.

Tissue or cell extraction samples prepared by chemical lysis buffer may cause unexpected ELISA results due to the impacts of certain chemicals.

7. Considering the possibility of mismatch between antigen from other resources and antibody in our kits (e.g., antibody targets conformational epitope rather than linear epitope), some native or recombinant proteins from other manufacturers may not be recognized by our products.

8. Due to factors including cell viability, cell number and sampling time, samples from cell culture supernatant may not be detected by the kit.

Recommend to use fresh samples for the test. If you store samples for long time, protein degradation and denaturalization may occur in those samples and finally lead to wrong results.

## Reagents Preparation

- Kindly use graduated containers to prepare the reagent. Please don't prepare the reagent directly in the Diluent vials provided in the kit.
- Bring all reagents to room temperature (18-25°C) before use for 30 mins.
- Prepare fresh standard for each assay. Use within 4 hours and discard after use
- It is not allowed to make serial dilution in the wells directly.
- Distilled Water is recommended. Contaminated water or container for reagents preparation will affect the test result.

**Biotin-conjugate (1x)** - Centrifuge the vial before opening.

Biotin-conjugate requires a 100-fold dilution. A suggested 100-fold dilution is 10 µl of Biotin-conjugate + 990 µl of Biotin-conjugate Diluent.

HRP-avidin (1x) - Centrifuge the vial before opening.

HRP-avidin requires a 100-fold dilution. A suggested 100-fold dilution is 10 μl of HRP-avidin + 990 μl of HRP-avidin Diluent.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Wash Buffer (1x)- If crystals have formed in the concentrate, warm them up

to room temperature and mix them gently until they get completely dissolved.

Dilute 20 ml of Wash Buffer Concentrate (25 x) into deionized or distilled

water to prepare 500 ml of Wash Buffer (1 x).

**BIOMATIK** 

Standard

Centrifuge the standard vial at 6000-10000rpm for 30s.

Dilute the **Standard(10x)** with **Sample Diluent**. A suggested 10-fold dilution

is 30 µl of Standard(10x) + 270 µl of Sample Diluent. This diluted Standard

(S7) serves as the high standard (80 ng/ml).

Do not substitute other diluents. Mix the standard to ensure complete dilution

and allow the standard to sit for a minimum of 15 mins with gentle agitation

prior to making dilutions.

Pipette 300 µl Sample Diluent into each tube (S0-S6). Use the diluted

solution to produce a 4-fold dilution series (below). Mix each tube thoroughly


before the next transfer. Sample Diluent serves as the zero standard (0

na/ml).



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140. (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

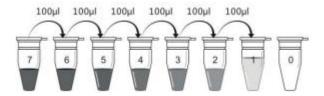

Tube	<b>S7</b>	<b>S</b> 6	S5	S4	S3	S2	S1	SO
ng/ml	80	20	5	1.25	0.312	0.078	0.019	0

#### **Key Notes**

- Do not mix or substitute reagents with those from other lots or sources
- If samples generate higher values than the highest standard, dilute the samples with Sample Diluent and repeat the assay.
- Any variation in Sample Diluent, operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding.
- This assay is designed to eliminate interference by soluble receptors, binding proteins, and other factors present in biological samples. Until all factors have been tested in the Immunoassay, the possibility of interference cannot be excluded.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Precautions

**BIOMATIK** 

The Stop Solution provided with this kit is an acid solution. Wear eve. hand, face.

and clothing protection when using this material.

Assay Procedures

Bring all reagents and samples to room temperature before use.

Centrifuge the sample again after thawing before the assay. It is

recommended to assay all samples and standards in duplicate.

1. Prepare all reagents, working standards, and samples as directed in

the previous sections.

2. Determine the number of wells to be used and put remaining wells and the

desiccant back into the pouch and seal the ziploc, store unused wells at 4°C.

3. Set a Blank well without any solution.

4. Add 50µl Standard or Sample per well.

5. Add 50µl Biotin-conjugate (1x) to each well(not to Blank well). Seal plate

with a new adhesive strip. Incubate at 37°C for 40 mins. (Biotin-conjugate

(1x) may appear cloudy. Warm up to room temperature and mix gently until

solution appears uniform.)

6. Aspirate each well and wash, repeating the process two times for a total of

three washes. Wash by filling each well with Wash Buffer (200µl) using a

squirt bottle, multi-channel pipette, manifold dispenser, or autowasher, and

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

let it stand for 2 mins, complete removal of liquid at each step is essential to

good performance. After the last wash, remove any remaining Wash Buffer

by aspirating or decanting. Invert the plate and blot it against clean paper

towels.

7. Add 100µl **HRP-avidin (1x)** to each well(not to Blank well). Seal plate with

a new adhesive strip. Incubate at 37°C for 40 mins.

8. Repeat the aspiration/wash process for five times as in step 6.

9. Add 90ul TMB Substrate to each well. Incubate at 37°C for 20 mins.

Protect from light.

**BIOMATIK** 

10. Add 50µl Stop Solution to each well, gently tap the plate to ensure

thorough mixing.

11. Determine the optical density of each well within 5 mins, using a

microplate reader set to 450 nm. If wavelength correction is available, set to

540 nm or 570 nm. Subtract readings at 540 nm or 570 nm from the readings

at 450 nm. This subtraction will correct for optical imperfections in the plate.

Readings made directly at 450 nm without correction may be higher and less

accurate

\*Samples may require dilution. Please refer to Sample Preparation.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Notes:

The final experimental results will be closely related to validity of 1

products, operation skills of end users and the experimental

environments.

**BIOMATIK** 

2. Samples or reagents addition: Please use the freshly prepared

Standard. Please carefully add samples to wells and mix gently to avoid

foaming. Do not touch the well wall as possible. For each step in the

procedure, total dispensing time for addition of reagents or samples to

the assay plate should not exceed 10 mins. This will ensure equal

elapsed time for each pipetting step, without interruption. Duplication of

all standards and specimens, although not required, is recommended. To

avoid cross-contamination, change pipette tips between additions of each

standard level, between sample additions, and between reagent additions.

Also, use separate reservoirs for each reagent.

3. Incubation: To ensure accurate results, proper adhesion of plate

sealers during incubation steps is necessary. Do not allow wells to sit

uncovered for extended periods between incubation steps. Once

reagents have been added to the well strips, DO NOT let the strips DRY

at any time during the assay. Incubation time and temperature must be

observed.





Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com



- Washing: The wash procedure is critical. Complete removal of 4. liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Solution by aspirating or decanting and remove any drop of water and fingerprint on the bottom of the plate. Insufficient washing will result in poor precision and falsely elevated absorbance reading. When using an automated plate washer, adding a 30 second soak period following the addition of wash buffer, and/or rotating the plate 180 degrees between wash steps may improve assay precision.
- 5. Controlling of reaction time: Observe the change of color after adding TMB Substrate (e.g. observation once every 10 mins), TMB Substrate should change from colorless or light blue to gradations of blue. If the color is too deep, add Stop Solution in advance to avoid excessively strong reaction which will result in inaccurate absorbance reading.
- 6. TMB Substrate is easily contaminated. TMB Substrate should remain colorless or light blue until added to the plate. Please protect it from light.
- 7 Stop Solution should be added to the plate in the same order as the TMB Substrate. The color developed in the wells will turn from blue to yellow upon addition of the Stop Solution. Wells that are green in color indicate that Stop Solution has not mixed thoroughly with TMB Substrate.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Calculation of Results

BIOMATIK

Using the professional soft "Curve Expert" to make a standard

curve is recommended, which can be downloaded from our web.

Average the duplicate readings for each standard and sample and

subtract the average optical density of Blank.

Create a standard curve by reducing the data using computer software

capable of generating a four parameter logistic (4-PL) curve-fit. As an

alternative, construct a standard curve by plotting the mean absorbance

for each standard on the x-axis against the concentration on the y-axis

and draw a best fit curve through the points on the graph. The data may

be linearized by plotting the log of the GHbA1c concentrations versus the

log of the O.D. and the best fit line can be determined by regression

analysis. This procedure will produce an adequate but less precise fit of

the data

If samples have been diluted, the concentration read from the standard

curve must be multiplied by the dilution factor.