



Email: info@biomatik.com https://www.biomatik.com

# **User Manual**

Catalog Number: EKC40095

Product Name: Human monocyte chemotactic protein 1/monocyte

chemotactic and activating factor, MCP-1/MCAF ELISA kit

Detection Range:6.25 pg/ml -400 pg/ml

Intended Use: For quantitative determination of Human monocyte chemotactic protein 1/monocyte chemotactic and activating factor (MCP-1/MCAF) concentrations in serum, plasma, urine, saliva, cerebrospinal fluid (CSF), ascitic fluid, cell culture supernates, tissue homogenates.

**Precautions:** For research use only. Not for use in diagnostic procedures.

Manual Version: 202301V1

# Storage:

Unopened kit	6 months when stored at 2 - 8°C.
Opened Kit	May be stored up to 1 month at 2 - 8°C. Keep it in sealed aluminum foil bag and avoid moisture.

The product manual may be updated as a result of continuous improvements.

Always refer to the hard copy manual included in the kit for your experiment.





Email: info@biomatik.com https://www.biomatik.com

### Kit Components

Reagents	Quantity		
Assay plate (12 x 8 coated Microwells)	1(96 wells)		
Standard (Freeze dried)	2		
Biotin-antibody (100 x concentrate)	1 x 120 µl		
HRP -avidin (100 x concentrate)	1 x 120 µl		
Biotin-antibody Diluent	1 x 15 ml		
HRP-avidin Diluent	1 x 15 ml		
Sample Diluent	1 x 50 ml		
Wash Buffer (25 x concentrate)	1 x 20 ml		
TMB Substrate	1 x 10 ml		
Stop Solution	1 x 10 ml		
Adhesive Strip (For 96 wells)	4		
Instruction manual	1		

#### **Working Principle**

This assay adopts quantitative sandwich enzyme immunoassay technique. Antibody specific to MCP-1/MCAF has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and then MCP-1/MCAF is bound by the immobilized antibody. After removing unbound substances, a biotin-conjugated antibody specific to MCP-1/MCAF is added to the wells. After washing, avidin conjugated Horseradish Peroxidase (HRP) is added to the wells. Following a wash to remove unbound avidin-enzyme reagent, then a substrate solution is added to the wells and later color develops in proportion to the amount of MCP-1/MCAF bound in the initial step. At last, measure the intensity of color after stopping color development.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Sensitivity

**BIOMATIK** 

The minimum detectable dose of mouse MCP-1/MCAF is less than 1.56

pg/ml. The sensitivity of this assay, or Lower Limit of Detection (LLD) was

defined as the lowest protein concentration that could be differentiated

from zero. It was determined by the mean O.D value of 20 replicates of

the zero standard added by their three standard deviations.

Specificity

This assay has high sensitivity and excellent specificity. No significant

cross-reactivity or interference between Human MCP-1/MCAF and

analogues was observed. Limited by current skills and knowledge, it is

impossible for us to complete the cross-reactivity detection between

Human MCP-1/MCAF and all the analogues, therefore, cross reaction

may still exist.

Precision

Intra-assay Precision (Precision within an assay): CV%<8%

Three samples of known concentration were tested twenty times on one

plate to assess.





Email: info@biomatik.com
https://www.biomatik.com

# Inter-assay Precision (Precision between assays): CV%<10%

Three samples of known concentration were tested in twenty assays to assess.

### Other Supplies Required

- Microplate reader capable of measuring absorbance at 450 nm,
   with the correction wavelength set at 540 nm or 570 nm.
- An incubator which can provide stable incubation conditions up to 37°C±0.5°C.
- Squirt bottle, manifold dispenser, or automated microplate washer.
- Absorbent paper for blotting the microtiter plate.
- 100ml and 500ml graduated cylinders.
- Deionized or distilled water.
- Pipettes and pipette tips.
- Test tubes for dilution.

Email: info@biomatik.com https://www.biomatik.com

**BIOMATIK** 

#### Sample Collection & Storage

• Serum Use a serum separator tube (SST) and allow samples to clot for two hours at room temperature or overnight at 4°C before centrifugation at 1000 ×g for 15 mins. Remove serum and assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles.

- Plasma Collect plasma using EDTA or heparin as an anticoagulant.
   Within 30 mins after collecting samples, centrifuge samples at 1000 x g, 2
   8°C, for 15 mins. Assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles.
- Tissue Homogenates Rinse 100mg tissue with 1X PBS, homogenize it in 1 ml of 1X PBS and store it overnight at -20°C. The cell membranes break after two freeze-thaw cycles, then centrifuge homogenates at 5000 x g, 2 8°C for 5 mins. Take the supernate for assay. Alternatively, aliquot and store the supernate at -20°C or -80°C. Centrifuge samples again after thawing. Avoid repeated freeze-thaw cycles.
- **Urine** Use a sterile container to collect urine samples. Remove any particulates by centrifugation for 15 minutes at 1000xg, 2 8°C and assay immediately or aliquot and store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles. Centrifuge again before assaying to remove

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

any additional precipitates that may appear after storage.

BIOMATIK

Saliva Remove particulates by centrifugation for 10 minutes at

4000x g at 2-8°C and assay immediately or aliquot and store samples at -

20°C. Avoid repeated freeze-thaw cycles. Recommend to use fresh

saliva samples.

Cerebrospinal Fluid (CSF) Remove particulates by centrifugation

for 10 minutes at 4000x g at 2-8°C and assay immediately or aliquot and

store samples at -20°C. Avoid repeated freeze-thaw cycles.

**Ascitic Fluid** Remove particulates by centrifugation for 10 minutes

at 4000x g at 2-8°C and assay immediately or aliquot and store samples

at -20°C. Avoid repeated freeze-thaw cycles.

Cell Culture Supernates Remove particulates by centrifugation for

15 minutes at 1000 x a. 2 - 8°C and assay immediately or aliquot and

store samples at -20°C or -80°C. Avoid repeated freeze-thaw cycles.

Notes:

1. Biomatik is only responsible for the kit itself, not for the samples

consumed during the assay. The user need to calculate the possible

amount of the samples to be used in the whole test. Please reserve

sufficient samples in advance.



Email: info@biomatik.com
https://www.biomatik.com

**BIOMATIK** 

2. Samples to be used within 5 days may be stored at 2-8°C, otherwise, samples must be stored at -20°C (≤1month) or -80°C (≤2month) to avoid contamination and loss of bioactivity.

3. Grossly hemolyzed samples are not suitable to use.

4. It would be necessary to run a preliminary experiment for validation, if the samples are not indicated in the manual.

5. Please predict the concentration before assaying. If results were not within the range of the standard curve, users would need to estimate the optimal sample dilutions for their particular experiments.

Tissue or cell extraction samples prepared by chemical lysis buffer may cause unexpected ELISA results due to the impacts of certain chemicals.

7. Considering the possibility of mismatch between antigen from other resources and antibody in our kits (e.g., antibody targets conformational epitope rather than linear epitope), some native or recombinant proteins from other manufacturers may not be recognized by our products.

Due to factors including cell viability, cell number and sampling time, samples from cell culture supernatant may not be detected by the kit.

Recommend to use fresh samples for the test. If you store samples for long time, protein degradation and denaturalization may occur in those samples and finally lead to wrong results.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Reagents Preparation

BIOMATIK

Kindly use graduated containers to prepare the reagent.

Please don't prepare the reagent directly in the Diluent vials in the

kit.

Bring all reagents to room temperature (18-25°C) before use for 30

mins.

Prepare fresh standard for each assay. Use it within 4 hours and

discard it after use

It is not allowed to make serial dilution in the wells directly.

Please carefully reconstitute Standards according to the instruction,

and avoid foaming and mix gently until the crystals have completely

dissolved. To minimize imprecision caused by pipetting, use small volumes

and ensure that pipettors are calibrated. It is recommended to suck more

than 10<sub>ul</sub> for once pipetting.

Distilled Water is recommended. Contaminated water or container for

reagents preparation will affect the test result.

**Biotin-antibody (1x)-** Centrifuge the vial before opening.

Biotin-antibody requires a 100-fold dilution. A suggested 100-fold dilution is

10 µl of Biotin-antibody+ 990 µl of Biotin-antibody Diluent.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

HRP-avidin (1x)- Centrifuge the vial before opening.

BIOMATIK

HRP-avidin requires a 100-fold dilution. A suggested 100-fold dilution is 10

μl of HRP-avidin + 990 μl of HRP-avidin Diluent.

Wash Buffer (1x)- If crystals have formed in the concentrate, warm up to

room temperature and mix gently until the crystals have completely dissolved.

Dilute 20 ml of Wash Buffer Concentrate (25 x) into deionized or distilled

water to prepare 500 ml of Wash Buffer (1 x).

Standard

Centrifuge the standard vial at 6000-10000 rpm for 30s.

Reconstitute the **Standard** with 1.0 ml of **Sample Diluent**. Do not substitute

other diluents. This reconstitution produces a stock solution of 400 pg/ml. Mix

the standard to ensure complete reconstitution and allow the standard to stay

for a minimum 15 mins with gentle agitation prior to making dilutions.

Pipette 250 ul Sample Diluent into each tube (S0-S6). Use the stock

solution to produce a 2-fold dilution series (below). Mix each tube thoroughly


before the next transfer. The undiluted Standard serves as the high standard

(400 pg/ml). **Sample Diluent** serves as the zero standard (0 pg/ml).





Email: info@biomatik.com https://www.biomatik.com

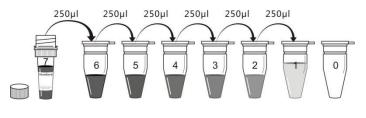

Tube	S7	S6	S5	S4	S3	S2	S1	S0
pg/ml	400	200	100	50	25	12.5	6.25	0

# **Key Notes**

- Do not mix or substitute reagents with those from other lots or sources.
- If samples generate higher values than the highest standard, dilute the samples with Sample Diluent and repeat the assay.
- Any variation in Sample Diluent, operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding.
- This assay is designed to eliminate interference by soluble receptors, binding proteins, and other factors present in biological samples. Until all factors have been tested in the Immunoassay, the possibility of interference cannot be excluded.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Precautions

The Stop Solution provided with this kit is an acid solution. Wear eve. hand, face.

and clothing protection when using this material.

BIOMATIK

Assay Procedures

Bring all reagents and samples to room temperature before use.

Centrifuge the sample again after thawing before the assay. It is

recommended to assay all samples and standards in duplicate.

1. Prepare all reagents, working standards, and samples as directed in

the previous sections.

2. Refer to the Assay Layout Sheet to determine the number of wells to be

used and put any remaining wells and the desiccant back into the pouch and

seal the ziploc, store unused wells at 4°C.

3. Add 100µl of standard and sample per well. Cover with the adhesive strip

provided. Incubate for 2 hours at 37°C. A plate layout is provided to record

standards and samples assayed.

4. Remove the liquid of each well, don't wash.

5. Add 100µl of Biotin-antibody (1x) to each well. Cover with a new

adhesive strip. Incubate for 1 hour at 37°C. (Biotin-antibody (1x) may

appear cloudy. Warm up to room temperature and mix gently until solution

appears uniform.)



Tel: (519) 489-7195, (800) 836-8089

Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

**BIOMATIK** 

6. Aspirate each well and wash, repeating the process two times for a total of

three washes. Wash by filling each well with Wash Buffer (200µl) using a

squirt bottle, multi-channel pipette, manifold dispenser, or autowasher, and

let it stand for 2 mins, complete removal of liquid at each step is essential to

good performance. After the last wash, remove any remaining wash Buffer

by aspirating or decanting. Invert the plate and blot it against clean paper

towels.

7. Add 100ul of **HRP-avidin (1x)** to each well. Cover the microtiter plate with

a new adhesive strip. Incubate for 1 hour at 37°C.

8. Repeat the aspiration/wash process for five times as in step 6.

Add 90ul of TMB Substrate to each well. Incubate for 15-30 mins at 37°C.

Protect from light.

10. Add 50µl of **Stop Solution** to each well, gently tap the plate to ensure

thorough mixing.

11. Determine the optical density of each well within 5 mins, using a

microplate reader set to 450 nm. If wavelength correction is available, set to

540 nm or 570 nm. Subtract readings at 540 nm or 570 nm from the readings

at 450 nm. This subtraction will correct for optical imperfections in the plate.

Readings made directly at 450 nm without correction may be higher and less

accurate.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

**BIOMATIK** 

Notes:

The final experimental results will be closely related to validity of 1

the products, operation skills of the end users and the experimental

environments.

2. Samples or reagents addition: Please use the freshly prepared

Standard. Please carefully add samples to wells and mix gently to avoid

foaming. Do not touch the well wall as possible. For each step in the

procedure, total dispensing time for addition of reagents or samples to

the assay plate should not exceed 10 mins. This will ensure equal

elapsed time for each pipetting step, without interruption. Duplication of

all standards and specimens, although not required, is recommended. To

avoid cross-contamination, change pipette tips between additions of each

standard level, between sample additions, and between reagent additions.

Also, use separate reservoirs for each reagent.

3. Incubation: To ensure accurate results, proper adhesion of plate

sealers during incubation steps is necessary. Do not allow wells to sit

uncovered for extended periods between incubation steps. Once

reagents have been added to the well strips, DO NOT let the strips DRY

at any time during the assay. Incubation time and temperature must be

observed.





Email: info@biomatik.com https://www.biomatik.com

4. Washing: The wash procedure is critical. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Solution by aspirating or decanting and remove any drop of water and fingerprint on the bottom of the plate. Insufficient washing will result in poor precision and falsely elevated absorbance reading. When using an automated plate washer, adding a 30 second soak period following the addition of wash buffer, and/or rotating the plate 180 degrees between wash steps may improve assay precision.

- 5. Controlling of reaction time: Observe the change of color after adding TMB Substrate (e.g. observation once every 10 mins), TMB Substrate should change from colorless or light blue to gradations of blue. If the color is too deep, add Stop Solution in advance to avoid excessively strong reaction which will result in inaccurate absorbance reading.
- TMB Substrate is easily contaminated. TMB Substrate should remain colorless or light blue until added to the plate. Please protect it from light.
- 7. Stop Solution should be added to the plate in the same order as the TMB Substrate. The color developed in the wells will turn from blue to yellow upon addition of the Stop Solution. Wells that are green in color indicate that Stop Solution has not mixed thoroughly with TMB Substrate.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Calculation of Results

BIOMATIK

Using the professional soft "Curve Expert" to make a standard

curve is recommended, which can be downloaded from our web.

Average the duplicate readings for each standard and sample and

subtract the average zero standard optical density.

Create a standard curve by reducing the data using computer software

capable of generating a four parameter logistic (4-PL) curve-fit. As an

alternative, construct a standard curve by plotting the mean absorbance

for each standard on the x-axis against the concentration on the y-axis

and draw a best fit curve through the points on the graph. The data may

be linearized by plotting the log of the MCP-1/MCAF concentrations

versus the log of the O.D. and the best fit line can be determined by

regression analysis. This procedure will produce an adequate but less

precise fit of the data.

If samples have been diluted, the concentration read from the standard

curve must be multiplied by the dilution factor.