

Email: info@biomatik.com
https://www.biomatik.com

# **User Manual**

Catalog Number: EKC40874

Product Name: Human placental protein 13(PP13) ELISA Kit

Detection Range: 3.9 pg/ml-1000 pg/ml

Intended Use: For quantitative determination of human placental protein 13

(PP13) concentrations in serum, urine.

Precautions: For research use only. Not for use in diagnostic procedures.

Manual Version: 202301V1

## Storage:

Unopened kit	6 months when stored at 2 - 8°C.
Opened Kit	May be stored up to 1 month at 2 - 8°C. Keep it in sealed aluminum foil bag and avoid moisture.

The product manual may be updated as a result of continuous improvements.

Always refer to the hard copy manual included in the kit for your experiment.





Email: info@biomatik.com https://www.biomatik.com

### Kit Components

Reagents	Quantity	
Assay Plate (12 x 8 coated Microwells)	1(96 wells)	
Standard	2 x 250 µl	
Antibody (100 x concentrate)	1 x 60 µl	
HRP-conjugate (100 x concentrate)	1 x 120 µl	
Antibody Diluent	1 x 10 ml	
HRP-conjugate Diluent	1 x 20 ml	
Sample Diluent	2 x 20 ml	
Wash Buffer (25 x concentrate)	1 x 20 ml	
TMB Substrate	1 x 10 ml	
Stop Solution	1 x 10 ml	
Adhesive Strip (For 96 wells)	4	
Instruction Manual	1	

## **Working Principle**

This assay employs the competitive inhibition enzyme immunoassay technique. The microtiter plate provided in this kit has been pre-coated with an antigen. Standards or samples are added to the appropriate microtiter plate wells with antibody specific for PP13 and Horseradish Peroxidase (HRP) conjugated goat-anti-rabbit antibody. The competitive inhibition reaction is launched between with pre-coated PP13 and PP13 in samples. A substrate solution is added to the wells and the color develops in opposite to the amount of PP13 in the sample. At last, measure the intensity of color after stopping color development.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Sensitivity

**BIOMATIK** 

The minimum detectable dose of human PP13 is less than 3.9 pg/ml. The

sensitivity of this assay, or Lower Limit of Detection (LLD) was defined as

the lowest human PP13 concentration that could be differentiated from

zero.

Specificity

This assay has high sensitivity and excellent specificity. No significant

cross-reactivity or interference between human PP13 and analogues was

observed. Limited by current skills and knowledge, it is impossible for us

to complete the cross-reactivity detection between human PP13 and all

the analogues, therefore, cross reaction may still exist.

Precision

Intra-assay Precision (Precision within an assay): CV%<12%

Three samples of known concentration were tested twenty times on one

plate to assess.

Inter-assay Precision (Precision between assays): CV%<15%

Three samples of known concentration were tested in twenty assays to

assess.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

**BIOMATIK** 

Other Supplies Required

Microplate reader capable of measuring absorbance at 450 nm,

with the correction wavelength set at 600 nm or 630 nm.

An incubator which can provide stable incubation conditions up to

37°C±0.5°C.

Squirt bottle, manifold dispenser, or automated microplate washer.

Absorbent paper for blotting the microtiter plate.

100ml and 500ml graduated cylinders.

Deionized or distilled water.

Pipettes and pipette tips.

Test tubes for dilution.

Sample Collection & Storage

**Serum** Use a serum separator tube (SST) and allow samples to

clot for two hours at room temperature or overnight at 4°C before

centrifugation at 1000 ×g for 15 mins. Remove serum and assay

immediately or aliquot and store samples at -20°C or -80°C. Avoid

repeated freeze-thaw cycles.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Urine Use a sterile container to collect urine samples. Remove any

particulates by centrifugation for 15 minutes at 1000x g, 2 - 8°C and

assay immediately or aliquot and store samples at -20°C or -80°C. Avoid

repeated freeze-thaw cycles. Centrifuge again before assaying to remove

any additional precipitates that may appear after storage.

Notes:

1. Biomatik is only responsible for the kit itself, not for the samples

consumed during the assay. The user need to calculate the possible

amount of the samples to be used in the whole test. Please reserve

sufficient samples in advance.

**BIOMATIK** 

2. Samples to be used within 5 days may be stored at 2-8°C, otherwise,

samples must be stored at -20°C (≤1month) or -80°C (≤2month) to avoid

contamination and loss of bioactivity.

3. Grossly hemolyzed samples are not suitable to use.

4. It would be necessary to run a preliminary experiment for validation, if

the samples are not indicated in the manual.

5. Please predict the concentration before assaying. If results were not

within the range of the standard curve, users would need to estimate the

optimal sample dilutions for their particular experiments.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

**BIOMATIK** 

6. Tissue or cell extraction samples prepared by chemical lysis buffer

may cause unexpected ELISA results due to the impacts of certain

chemicals

7. Considering the possibility of mismatch between antigen from other

resources and antibody in our kits (e.g., antibody targets conformational

epitope rather than linear epitope), some native or recombinant proteins

from other manufacturers may not be recognized by our products.

8. Due to factors including cell viability, cell number and sampling time.

samples from cell culture supernatant may not be detected by the kit.

9. Recommend to use fresh samples for the test. If you store samples for

long time, protein degradation and denaturalization may occur in those

samples and finally lead to wrong results.

Reagents Preparation

Kindly use graduated containers to prepare the reagent.

Please don't prepare the reagent directly in the Diluent vials

provided in the kit.

Bring all reagents to room temperature (18-25°C) before use for 30

mins

Prepare fresh standard for each assay. Use within 4 hours and discard

after use.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

Making serial dilution in the wells directly is not permitted.

Distilled Water is recommended. Contaminated water or container for

reagents preparation will affect the test result.

**Antibody (1x)** - Centrifuge the vial before opening.

Antibody requires a 100-fold dilution. A suggested 100-fold dilution is 10 µl

of Antibody + 990 µl of Antibody Diluent.

**BIOMATIK** 

**HRP-conjugate (1x)** - Centrifuge the vial before opening.

HRP-conjugate requires a 100-fold dilution. A suggested 100-fold dilution is

10 µl of HRP-conjugate + 990 µl of HRP-conjugate Diluent.

Wash Buffer (1x)- If crystals have formed in the concentrate, warm them up

to room temperature and mix them gently until they get completely dissolved.

Dilute 20 ml of Wash Buffer Concentrate (25 x) into deionized or distilled

water to prepare 500 ml of Wash Buffer (1 x).

Standard

Centrifuge the standard vial at 6000-10000rpm for 30s before opening.


Allow the stock standard solution of 1000 pg/ml to sit for a minimum of 30

mins to warm to room temperature prior to making dilutions.



Email: info@biomatik.com https://www.biomatik.com

Pipette 150 µl **Sample Diluent** into each tube (S0-S4). Use the stock solution to produce a 4-fold dilution series (below). Mix each tube thoroughly before the next transfer. The undiluted Standard serves as the high standard (1000 pg/ml). **Sample Diluent** serves as the zero standard (0 pg/ml).

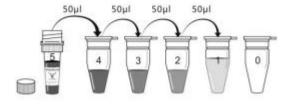

Tube	S5	S4	S3	S2	S1	S0
pg/ml	1000	250	62.5	15.6	3.9	0

# **Key Notes**

- Do not mix or substitute reagents with those from other lots or sources.
- If samples generate higher values than the highest standard, dilute the samples with Sample Diluent and repeat the assay.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

**BIOMATIK** 

Any variation in Sample Diluent, operator, pipetting technique,

washing technique, incubation time or temperature, and kit age can

cause variation in binding.

This assay is designed to eliminate interference by soluble

receptors, binding proteins, and other factors present in biological

samples. Until all factors have been tested in the Immunoassay, the

possibility of interference cannot be excluded.

Precautions

The Stop Solution provided with this kit is an acid solution. Wear eye, hand, face,

and clothing protection when using this material.

Assav Procedures

Bring all reagents and samples to room temperature before use.

Centrifuge the sample again after thawing before the assay. It is

recommended to assay all samples and standards in duplicate.

1. Prepare all reagents, working standards, and samples as directed in

the previous sections.

2. Determine the number of wells to be used and put remaining wells and the

desiccant back into the pouch and seal the ziploc, store unused wells at 4°C.

3. Set a **Blank** well without any solution.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

**BIOMATIK** 

4. Add 50µl Standard or Sample per well. Add 50µl Antibody (1x) to each well immediately (not to Blank well). Mix well with the pipette or shake the plate gently for 60 seconds. A plate layout is provided to record standards and samples assayed.

5. Seal plate with adhesive strip. Incubate at 37°C for 30 mins.

6. Aspirate each well and wash, repeating the process two times for a total of three washes. Wash by filling each well with Wash Buffer (250µl) using a squirt bottle, multi-channel pipette, manifold dispenser, or autowasher, and let it stand for 2 mins, complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.

7. Add 100ul HRP-conjugate (1x) to each well immediately (not to Blank well). Seal plate with adhesive strip. Incubate at 37°C for 30 mins.

8. Repeat the aspiration/wash process for five times as in step 6.

9. Add 90µl TMB Substrate to each well. Incubate at 37°C for 20 mins. Protect from light.

10. Add 50ul Stop Solution to each well, gently tap the plate to ensure thorough mixing.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

11. Determine the optical density of each well within 5 mins, using a

microplate reader set to 450 nm. If wavelength correction is available, set to

540 nm or 570 nm. Subtract readings at 540 nm or 570 nm from the readings

at 450 nm. This subtraction will correct for optical imperfections in the plate.

Readings made directly at 450 nm without correction may be higher and less

accurate.

**BIOMATIK** 

Notes:

The final experimental results will be closely related to validity of 1.

products, operation skills of end users and the experimental

environments.

2. Samples or reagents addition: Please use the freshly prepared

Standard, Please carefully add samples to wells and mix gently to avoid

foaming. Do not touch the well wall as possible. For each step in the

procedure, total dispensing time for addition of reagents or samples to

the assay plate should not exceed 10 mins. This will ensure equal

elapsed time for each pipetting step, without interruption. Duplication of

all standards and specimens, although not required, is recommended. To

avoid cross-contamination, change pipette tips between additions of each

standard level, between sample additions, and between reagent additions.

Also, use separate reservoirs for each reagent.





Email: info@biomatik.com



- 3. Incubation: To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary. Do not allow wells to sit uncovered for extended periods between incubation steps. Once reagents have been added to the well strips, DO NOT let the strips DRY at any time during the assay. Incubation time and temperature must be observed.
- 4. Washing: The wash procedure is critical. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Solution by aspirating or decanting and remove any drop of water and fingerprint on the bottom of the plate. Insufficient washing will result in poor precision and falsely elevated absorbance reading. When using an automated plate washer, adding a 30 second soak period following the addition of wash buffer, and/or rotating the plate 180 degrees between wash steps may improve assay precision.
- 5. Controlling of reaction time: Observe the change of color after adding TMB Substrate (e.g. observation once every 10 mins), TMB Substrate should change from colorless or light blue to gradations of blue. If the color is too deep, add Stop Solution in advance to avoid excessively strong reaction which will result in inaccurate absorbance reading.

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com

https://www.biomatik.com

TMB Substrate is easily contaminated. TMB Substrate should 6.

remain colorless or light blue until added to the plate. Please protect it

from light.

7. Stop Solution should be added to the plate in the same order as

the TMB Substrate. The color developed in the wells will turn from blue to

yellow upon addition of the Stop Solution. Wells that are green in color

indicate that Stop Solution has not mixed thoroughly with TMB Substrate.

Calculation of Results

**BIOMATIK** 

Using the professional soft "Curve Expert" to make a standard

curve is recommended, which can be downloaded from our web.

Average the duplicate readings for each standard and sample and

subtract the average optical density of Blank.

Create a standard curve by reducing the data using computer software

capable of generating a four parameter logistic (4-PL) curve-fit. As an

alternative, construct a standard curve by plotting the mean absorbance

for each standard on the x-axis against the concentration on the y-axis

and draw a best fit curve through the points on the graph. The data may

be linearized by plotting the log of the PP13 concentrations versus the log



Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515

Email: info@biomatik.com https://www.biomatik.com

of the O.D. and the best fit line can be determined by regression analysis.

This procedure will produce an adequate but less precise fit of the data.

If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.