

Biomatik Tel:(519) 489-7195,(800) 836-8089 Fax:(519) 231-0140,(877) 221-3515 Email:info@biomatik.com http://www.biomatik.com

RPU40107

Recombinant Adenylyl Cyclase Associated Protein 1 (CAP1)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)



Biomatik Tel:(519) 489-7195,(800) 836-8089 Fax:(519) 231-0140,(877) 221-3515 Email:info@biomatik.com http://www.biomatik.com

[PROPERTIES]

Source: Prokaryotic expression

Host: E.coli

Residues: Ala38~Lys149

Tags: N-terminal His Tag

Subcellular Location: Membrane

Purity: > 95%

Traits: Freeze-dried powder

Buffer formulation: 20mM Tris, 150mM NaCl, pH8.0, containing 0.01% SKL, 5% Trehalose.

Original Concentration: 500µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.8

Predicted Molecular Mass: 13.4kDa

Accurate Molecular Mass: 14kDa as determined by SDS-PAGE reducing conditions.

[<u>USAGE</u>]

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[<u>SEQUENCE</u>]



Biomatik Tel:(519) 489-7195,(800) 836-8089 Fax:(519) 231-0140,(877) 221-3515 Email:info@biomatik.com http://www.biomatik.com

AGA APYVQAFDSL LAGPVAEYLK ISKEIGGDVQ KHAEMVHTGL KLERALLVTA SQCQQPAENK LSDLLAPISE QIKEVITFRE KNRGSKLFNH LSAVSESIQA LGWVAMAPK

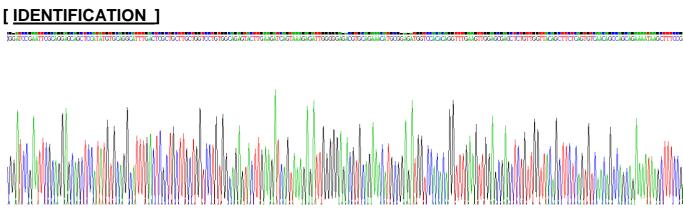


Figure . Gene Sequencing (extract)

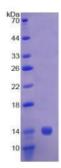


Figure. SDS-PAGE

[<u>IMPORTANT NOTE</u>]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.