

Biomatik Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515 Email: info@biomatik.com http://www.biomatik.com

RPU41008 50µg

**Recombinant Coronin 1A (CORO1A)** 

Organism Species: Homo sapiens (Human)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)



Biomatik Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515 Email: info@biomatik.com http://www.biomatik.com

# [ PROPERTIES ]

Source: Prokaryotic expression

Host: E.coli

Residues: Arg7~Glu204

Tags: N-terminal His Tag

Subcellular Location: Cytoplasm

**Purity:** > 95%

Traits: Freeze-dried powder

Buffer formulation: 100mMNaHCO<sub>3</sub>, 500mMNaCl, pH8.3, containing 0.01% SKL, 5%

Trehalose.

**Original Concentration:** 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.3

Predicted Molecular Mass: 25.1kDa

Accurate Molecular Mass: 25kDa as determined by SDS-PAGE reducing conditions.

# [<u>USAGE</u>]

Reconstitute in 100mM NaHCO3, 500mM NaCl (pH8.3) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

# [ STORAGE AND STABILITY ]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.




Biomatik Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515 Email: info@biomatik.com http://www.biomatik.com

### [SEQUENCE]

RSSK FRHVFGQPAK ADQCYEDVRV SQTTWDSGFC AVNPKFVALI CEASGGGAFL VLPLGKTGRV DKNAPTVCGH TAPVLDIAWC PHNDNVIASG SEDCTVMVWE IPDGGLMLPL REPVVTLEGH TKRVGIVAWH TTAQNVLLSA GCDNVIMVWD VGTGAAMLTL GPEVHPDTIY SVDWSRDGGL ICTSCRDKRV RIIE

# [ IDENTIFICATION ]

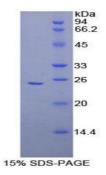

# [<u>IMPORTANT NOTE</u>]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.