

Biomatik

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515 Email: info@biomatik.com http://www.biomatik.com

RPU41409 100µg

**Recombinant Fibrinogen Alpha Chain (FGA)** 

Organism Species: Homo sapiens (Human)

Instruction manual

FOR RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)







### [PROPERTIES]

Source: Prokaryotic expression

Host: E.coli

Residues: Asp124~Pro214

Tags: N-terminal His Tag

Subcellular Location: Secreted

**Purity:** > 97%

Traits: Freeze-dried powder

Buffer formulation: 100mM NaHCO<sub>3</sub>, 500mM NaCl, pH8.3, containing 0.01% sarcosyl,

5%Trehalose.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.4

Predicted Molecular Mass: 17.4kDa

Accurate Molecular Mass: 14kDa as determined by SDS-PAGE reducing conditions.

#### [USAGE]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-1.0 mg/mL. Do not vortex.

#### [STORAGE AND STABILITY]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

#### [SEQUENCE]



**Biomatik** 

Tel: (519) 489-7195, (800) 836-8089 Fax: (519) 231-0140, (877) 221-3515 Email: info@biomatik.com http://www.biomatik.com

DNTYNRV SEDLRSRIEV LKRKVIEKVQ HIQLLQKNVR AQLVDMKRLE VDIDIKIRSC RGSCSRALAR EVDLKDYEDQ QKQLEQVIAK DLLP

## [ IDENTIFICATION ]

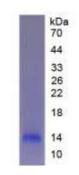


Figure. SDS-PAGE

# [ IMPORTANT NOTE ]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.