



**RPU42126 10µg**

**Recombinant Interferon Alpha/Beta Receptor 1 (IFNα/bR1)**

**Organism Species: *Homo sapiens* (Human)**

**Instruction manual**

**FOR IN VITRO USE AND RESEARCH USE ONLY  
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES**

12th Edition (Revised in Aug, 2016)



## [ PROPERTIES ]

**Source:** Prokaryotic expression

**Host:** *E.coli*

**Residues:** Lys28~Asn227

**Tags:** N-terminal His Tag

**Subcellular Location:** Membrane

**Purity:** > 90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 7.2

**Predicted Molecular Mass:** 24.2kDa

**Accurate Molecular Mass:** 26kDa as determined by SDS-PAGE reducing conditions.

## [ USAGE ]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

## [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ SEQUENCE ]

Lys28 ~ Asn227

```
                KNL KSPQKVEVDI IDDNFILRWN  
RSDSVGNTV FSFDYQKTGM DNWIKLSGCQ NITSTKCNFS SLKLVYEEI  
KLRIRAEKEN TSSWYEVSF TPFKAQIGP PEVHLEAEDK AIVIHISPGT  
KDSVMWALDG LSFTYSLVIW KNSSGVEERI ENIYSRHKIY KLSPETTYCL  
KVKAALLTSW KIGVYSPVHC IKTTVEN
```

## [ IDENTIFICATION ]

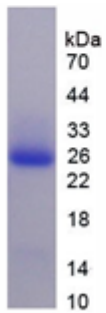


Figure. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.