



RPU43963 10µg

Recombinant Telomerase Reverse Transcriptase (TERT)

Organism Species: *Homo sapiens* (Human)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Arg787~Arg1084

Tags: N-terminal His and GST Tag

Subcellular Location: Nucleus, Chromosome

Purity: > 97%

Traits: Freeze-dried powder

Buffer formulation: 20mM Tris, 150mM NaCl, pH8.0, containing 1mM EDTA, 1mM DTT, 0.01%skl, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.8

Predicted Molecular Mass: 63.4kDa

Accurate Molecular Mass: 65kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

```

                                RDAV VIEQSSSLNE
ASSGLFDVFL RFMCHHAVRI RGKSYVQCQG IPQGSILSTL LCSLCYGDME
NKLFAGIRRD GLLLRLVDDF LLVTPHLTHA KTFLRTLVRG VPEYGCVVNL
RKTVVNFPVE DEALGGTAFV QMPAHGLFPW CGLLLDTRTL EVQSDYSSYA
RTSIRASLTF NRGFKAGRNM RRKLFGLVRL KCHSLFLDLQ VNSLQTVCTN
IYKILLQAY RFHACVLQLP FHQQVWKNPT FFLRVISDTA SLCYSILKAK
NAGMSLGAKG AAGPLPSEAV QWLCHQAFLL KLTR
  
```

[IDENTIFICATION]

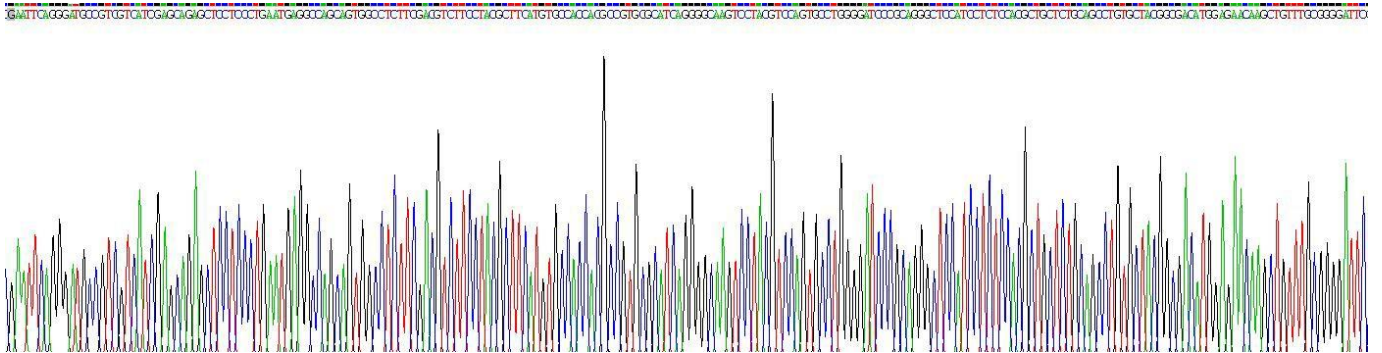


Figure . Gene Sequencing (extract)

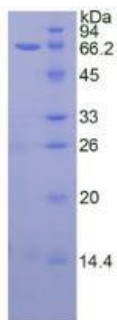


Figure. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.